


As6



(19)

Europäisches Patentamt  
European Patent Office  
Office européen des brevets

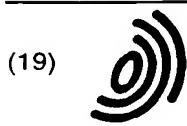

(11)

EP 0 877 092 A1

(12)

## EUROPEAN PATENT APPLICATION

(43) Date of publication:  
11.11.1998 Bulletin 1998/46

(51) Int Cl.<sup>6</sup>: C12Q 1/04

(21) Application number: 98302542.0

(22) Date of filing: 01.04.1998

(84) Designated Contracting States:  
AT BE CH CY DE DK ES FI FR GB GR IE IT LI LU  
MC NL PT SE  
Designated Extension States:  
AL LT LV MK RO SI

(30) Priority 08.04.1997 US 826844

(71) Applicant **Becton, Dickinson and Company**  
**Franklin Lakes, New Jersey 07417-1880 (US)**

(72) Inventors  
• **Strenkoski, Leon, F.**  
**Portland ME 04103-1613 (US)**

- **Schneider, Maureen, A.**  
**Farmington Hills, Michigan 49334 (US)**
- **Swamy, Shoba, C.**  
**Canton, Michigan 48188 (US)**
- **Nagar, Mandar, S.**  
**Ypsilanti Twp, Michigan 48197 (US)**

(74) Representative: **Austin, Hedley William**  
**Urquhart-Dykes & Lord**  
**Alexandra House**  
**1 Alexandra Road**  
**Swansea Wales SA1 5ED (GB)**

(54) **Rapid microorganism detection method**

(57) A rapid screening method for the detection and identification of target microorganisms in a sample which may contain both target microorganisms and competing non-target microorganisms includes the steps of enriching a sample in a mildly selective medium, incubating the sample in the growth medium for a predetermined amount of time, adding at least one inhibitor of the non-target microorganisms to the growth medium to discourage the growth of non-target micro-

organisms and to encourage the growth of target microorganisms. The method further includes the steps of incubating the sample in the growth medium including at least one inhibitor for a predetermined amount of time, performing biochemical assays specific for the identification of the target microorganisms, and detecting the presence of the target microorganisms in the sample. The invention further provides a system for the rapid detection of microorganisms in a sample.

**Description****TECHNICAL FIELD**

5 The present invention relates to a method for detecting the presence of microorganisms in a sample. More specifically, the present invention relates to a rapid method for detecting and identifying target microorganisms in a sample which may contain other competing microorganisms.

**BACKGROUND OF THE INVENTION**

10 Foods, pharmaceuticals, cosmetics, and water are routinely tested for microbial contamination with pathogenic microorganisms. A primary technique for screening for the presence of pathogens involves a series of media (nutrient fluid) transfers starting from a non-selective enrichment media and then transcending a portion of the primary enrichment medium to a selective media. This process allows for the initiation of growth of potentially injured microorganisms in the non-selective primary enrichment media and, once the microorganisms have been revived, a small quantity of the non-selective enrichment medium is transferred into the secondary (selective) media. This process, whose duration is often defined by human work patterns and the growth patterns of the microorganisms, can take several days to complete.

15 For example, culturing a food sample for the presence of *Salmonella* typically involves the addition of approximately 20-25 grams of sample, such as meat, into approximately 225 ml of primary enrichment media. The primary enrichment media is typically non-selective, such as Buffered Peptone Water (BPW) or Universal Pre-enrichment Broth (UPB), to allow for repair of injured microorganisms. The sample is then thoroughly mixed with the primary enrichment broth and incubated for 22-28 hours at approximately  $35^{\circ}\text{C} \pm 2^{\circ}\text{C}$ . Following this step, the sample is further selectively enriched in a growth promoting medium containing inhibitors that allow for the continued growth of a target organism, such as *Salmonella*, while simultaneously restricting the proliferation of most other competing microorganisms.

20 The most commonly used method for selective enrichment of the target organisms requires the transfer of approximately one milliliter of the primary enrichment media into two tubes containing ten milliliters of selective media such as Selenite Cystine Broth and Tetrathionate Broth (AOAC), respectively, and incubating these tubes for 22-28 hours at approximately  $35^{\circ}\text{C}$ . Other secondary enrichment media for *Salmonella* could include Rappaport-Vassiliadis Medium (RV) and Lauryl Tryptose Broth. This secondary enrichment step is generally followed by a detection step that can include plating the secondary enrichment broth onto solid media or by utilizing other more rapid methods, such as immunological assays or DNA probes (all media available from Difco Laboratories, Detroit, Michigan).

25 Many attempts have been made over the years to shorten the length of the primary enrichment step in order to reduce total assay time. Generally, it has been found that six (6) to twelve (12) hours of primary enrichment are required to obtain enough viable microorganisms for further testing. Some of these general methods are discussed below.

30 It has long been recognized that the potential for bacterial contamination exists in products which are designed for human consumption or use. Manufacturers and/or processors of these products destined for human consumption generally test these products to ensure the quality and safety of the products for human consumption. Such products can include raw meat, prepared food items, food preparation equipment, drinking water, bathing water, and other vectors in which microorganisms can reside which allows their contact or transfer to humans or animals.

35 Typically, organisms such as enteric bacteria, i.e., Enterobacteriaceae, such as *Salmonella* and *Escherichia coli*, and Gram positive organisms, such as *Staphylococcus* and *Enterococcus*, are the organisms tested for in products destined for human and animal consumption or use.

40 These microorganisms are generally present in the colon, intestines, or fecal matter of humans or animals. When food products such as poultry, red meat, seafood, eggs, or any foods which contain these products come into contact with fecal matter during handling or processing, the potential exists for contamination and subsequent transfer of these organisms to the end users or consumers, i.e., humans.

45 Not only can the presence of sufficient numbers of microorganisms cause the deterioration of a food product such as by causing spoilage, additionally, if consumed by a human or animal, can also cause disease. In the United States alone, the number of cases of food poisoning associated with the consumption of contaminated food products is conservatively estimated to be in the multi-millions per year. While most human cases of bacterial food poisoning only result in acute symptomatic disease which includes nausea, vomiting, diarrhea, chills, fever, and exhaustion; for those individuals such as infants, the elderly, pregnant women, neonates, and those with immunocompromised systems, death can occur.

50 The total economic loss attributable to bacterial food poisoning has been estimated to reach into the hundreds of millions of dollars each year due to lost productivity, increased use of the medical insurance system, and increased use of the medical provider system.

55 In order to prevent the transmission of foodborne bacterial pathogens, the manufacturers and/or processors of

food products routinely test samples of their products in order to identify contaminated products before the product is placed into the stream of commerce leading to human consumption. Because pathogenic microorganisms can be present in very small numbers in a food product which may contain a large number or variety of other pathogenic and/or non-pathogenic microorganisms, methods for the recovery and detection of sub-lethally injured pathogenic microorganisms have been developed. Pathogenic microorganisms are typically injured during the processing of the product destined for human consumption. That is, the pathogenic microorganisms may have undergone heating, freezing, contact with chemical additives, or mechanical processing steps which injure or debilitate the pathogenic microorganism present in the product.

A traditional method for recovering pathogenic microorganisms from a food sample involves five basic steps as set forth in United States Patent Number 5,145,786 to Bailey et al. The first step is a pre-enrichment wherein the food sample is enriched in a non-selective medium to restore injured bacterial cells to a stable physiological condition. The second step involves a selective enrichment wherein selective inhibitory reagents are added to the growth promoting medium to promote the growth of selected pathogenic microorganisms while restricting the proliferation of most other bacteria. The third step involves selectively plating a sample of the enriched growth medium onto a solid selective media to physically isolate pure, discrete colonies of the suspected bacterial pathogen. The next step involves biochemical testing or screening of the pure cultures obtained from the selective plating step in order to eliminate non-target or competing organisms. The last step involves a serological analysis of the pure culture of the suspected pathogenic microorganisms in order to specifically identify the pathogenic microorganism present in the sample. The major disadvantage of this conventional method is that it is very labor intensive and time consuming. The conventional method takes on the order of three to five days to obtain a positive or negative outcome for the presence or absence of the pathogenic microorganism in the analyzed sample and is both labor and resource intensive. This lengthy analysis is often not suitable for producers of perishable products because the producer must keep the product from the consumer until the test period is over and a negative test result is obtained.

In order to eliminate the lengthy time involved with the conventional method of detection, other methods have been developed which significantly reduce the time necessary to obtain a result. The REVEAL *Salmonella* Test System (Neogen Corp., Lansing, Michigan, U.S.P.N. 5,296,370) shortens the time frame necessary to obtain a negative result to approximately two days. On day one, a sample to be tested is inoculated into a non-inhibitory pre-enrichment media. After approximately two to four hours of incubation, highly selective inhibitors are added to selectively enrich for specific bacterial pathogens while inhibiting the proliferation of other microorganisms present in the sample. On day two, all of the samples are assayed by immunoassay in order to detect positive and negative samples. All positive samples are then tested further in order to confirm their identity. This method has the disadvantage that all of the samples must be analyzed by immunoassay in order to detect positive or negative cultures.

Other methods which are described and shown in Table 1, illustrate a number of commercially available test systems for the identification of *Salmonella* in a sample. All of the methods described in Table 1 require either the plating of selectively grown microorganisms onto a solid selective media to obtain a pure culture or require an immunological or nucleic acid based assay to be performed on every sample in order to obtain a positive or negative result for the presence or absence of *Salmonella* in the sample.

None of the methods described above or in Table 1 allows for a liquid enrichment of a sample to be tested for the presence of a microorganism, such as *Salmonella*, by adjusting the levels of *Salmonella* in the enrichment medium to be at or above the levels of contaminants (other microorganisms present in the sample) over a wide range of contaminant concentrations. Additionally, none of the above-described methods allow for liquid biochemical identification tests which have been fine-tuned to be reliably performed in mixed cultures. The above-described methods either require plating of an enriched culture onto a solid selective medium to obtain pure cultures prior to the performance of biochemical identification tests or require that all isolated samples be subjected to an immunoassay in order to identify the presence of a pathogenic microorganism such as *Salmonella*.

Therefore, it would be advantageous to have a rapid screening method for the detection and identification of target microorganisms, such as *Salmonella*, in a sample which may contain both target microorganisms and competing non-target microorganisms by utilizing both a liquid enrichment of sample to be tested for the presence of microorganism, such as *Salmonella*, by adjusting the levels of *Salmonella* in the enrichment medium to be at or above the levels of contaminants (other microorganisms present in the sample) over a wide range of contaminant concentrations and biochemical identification tests performed on potentially mixed cultures. It would further be advantageous to have a rapid method for the detection and identification of target microorganisms in a sample which may contain both target microorganisms and competing non-target microorganisms in which by utilizing biochemical tests, one would be able to rapidly screen out a large proportion of negative samples from further testing and thereby allow a manufacturer or processor of a product to release those products for consumption which have tested negative for the presence of pathogenic microorganisms.

The present invention not only provides a streamlined method for the identification and detection of target microorganisms in a sample, it yields the added benefits of allowing biochemical assays to be performed on potentially mixed

BEST AVAILABLE COPY

cultures and also provides a method to rapidly screen out a large percentage of negative samples from further testing. The present invention provides a method which can be adapted for use with a number of rapidly growing Gram negative (-) organisms as well as rapidly growing Gram positive (+) microorganisms.

## SUMMARY OF THE INVENTION

A method according to the present invention for the detection and identification of target microorganisms in a sample which may contain both target microorganisms and competing non-target microorganisms includes the steps of enriching a sample in a mildly selective medium, incubating the sample in the growth medium for a predetermined amount of time, adding at least one inhibitor of the non-target microorganisms to the growth medium to discourage the growth of non-target microorganisms and to encourage the growth of target microorganisms, incubating the sample in the growth medium for a predetermined amount of time, performing biochemical assays specific for the identification of the target microorganisms, and detecting the presence of the target microorganism in the sample. A system for the detection of microorganisms is also disclosed.

## BRIEF DESCRIPTION OF THE DRAWINGS

Other advantages of the present invention will be readily appreciated as the same becomes better understood by reference to the following detailed description when considered in connection with the accompanying drawings wherein:

Figure 1 is a graph illustrating the growth of heat injured *S. typhimurium* in various pre-enrichment media;

Figure 2 is a graph illustrating the growth over time of various bacterial samples in modified-TSB media with no inhibitors added

Figure 3 illustrates the growth over time of various bacterial samples in a lactose broth control medium;

Figure 4 is a graph illustrating the growth over time of various bacterial samples in modified-TSB media with the inhibitor combination brilliant green/deoxycholate;

Figure 5 is a graph illustrating the growth over time of various bacterial samples in lactose broth with the inhibitor combination brilliant green/deoxycholate;

Figure 6 is a graph illustrating the growth over time of various bacterial samples in modified-TSB with the inhibitor combination brilliant green/novobiocin;

Figure 7 is a graph illustrating the growth over time of various bacterial specimens in modified-TSB media with the inhibitor tergitol 4;

Figure 8 is a graph illustrating the growth over time of various bacterial specimens in lactose broth with the growth inhibitor tergitol 4;

Figure 9 is a graph illustrating the growth over time of various bacterial specimens in modified-TSB with the inhibitor combination tergitol 4/brilliant green;

Figure 10 is a graph illustrating the growth curve for *S. typhimurium* in modified-TSB media with novobiocin wherein the concentration of novobiocin was varied in order to determine its effect on growth; and

Figure 11 is a graph illustrating the growth curve for *S. brandenberg* in modified-TSB media with novobiocin wherein the concentration of novobiocin was varied in order to determine its effect on the growth of the *S. brandenberg*.

## DETAILED DESCRIPTION OF THE INVENTION

Generally, the present invention provides a rapid screening method for the detection and identification of rapidly growing target microorganisms in a sample which may contain both target microorganisms and competing non-target microorganisms. The method includes the steps of enriching a sample in a growth medium, incubating the sample in the growth medium for a predetermined amount of time, adding at least inhibitor of the non-target microorganisms to the growth medium to discourage the growth of non-target microorganisms and to encourage the growth of the target

microorganisms, and incubating the sample in the medium which includes an inhibitor(s) for a predetermined amount of time. The present invention also includes the steps of performing biochemical assays specific for the identification of the target microorganisms and then detecting the presence of the target microorganisms in the sample.

Target organisms are those organisms which are being specifically tested for and/or selected for using the method of the present invention. For example, in testing a food sample, the rapidly growing target microorganisms which could be tested for could include, but are not limited to, *Salmonella* and/or *E. coli*. Since these are particular pathogens which are often found in foods, they are often "targeted" for testing in order to determine if they are present in a food sample.

In general, the rapidly growing target microorganisms have a generation time of less than one hour. Preferably, the rapidly growing target microorganisms which can be rapidly detected by the method of the present invention will have a generation time of less than thirty minutes such as *Salmonella* and *E. coli*. Generation times for these organisms can range from fifteen to forty minutes depending on the type of media in which the organisms are grown. A rich media, such as TSB, can support lower (shorter) generation times whereas a weak media, such as minimal media, can or may only be able to support higher (longer) generation times.

The target microorganisms which can be detected and identified using the method of the present invention can include any rapidly growing Gram (-) organisms such as the Enterobacteriaceae including *E. coli*, and rapidly growing Gram (+) organisms such as *Staphylococcus aureus* and *Enterococcus faecalis*. The target microorganism can be detected and identified by tailoring the specific inhibitors and biochemical reagents to the target organism(s) to allow for the use of these inhibitors and biochemical tests in mixed culture as discussed below in greater detail.

In order to screen a sample for the presence of the target microorganism in the sample, a suitable sample i.e., a food sample, pharmaceutical sample, or other, is obtained. A sample of suitable size and weight is placed in a suitable container containing a suitable pre-enrichment media.

In a preferred embodiment of the present invention, the container is constructed in the form of bag. The bag can be constructed of any suitable resilient plastic material such as polyester/polyethylene laminate. Because of the nature of the samples which may be tested, such as a food sample, the container should be constructed of a resilient material so that it can withstand homogenizing or "stomaching" of the sample disposed within the bag in order to expose the microorganisms located within the sample. The container can also be made of any other suitable material known to those skilled in the art.

The container can also include a closure such as a "zip-lock" type closure to prevent spillage of the contents of the container and, further, to prevent contamination of the contents of the container. The container can also include a filter which allows fluid to pass therethrough while preventing particular matter, such as food particles, from passing therethrough.

Following dispensing of the sample in the container of pre-enrichment media, the sample, can be stomached or homogenized within the bag utilizing means well known to those skilled in the art. (Stomacher Lab-Blender, Seward Medical, London, U.K.). The sample is stomached within the bag containing the pre-enrichment media for a period of approximately thirty seconds to five minutes. The container is then incubated at a temperature suitable for growth of the target microorganism. For example, if the target organism is *Salmonella*, then the container could be incubated at approximately  $35^{\circ}\text{C} \pm 2^{\circ}\text{C}$ . The incubation in the pre-enrichment media (pre-enrichment step) is carried out to resuscitate sub-lethally injured microorganisms that are present in the sample.

The pre-enrichment step may be carried out for a period of approximately two to six hours in duration. Preferably, the pre-enrichment step is carried out for period of approximately four hours in duration.

The pre-enrichment media can include any general, non-selective and/or mildly selective growth medium known to those skilled in the art. This can include media such as REVIVE ((Neogen Corp., Lansing, Michigan) TSB media (Difco, Detroit, Michigan) or other similar type media. The pre-enrichment medium can be any non-selective liquid media to a mildly selective growth medium allowing for the rapid recovery and growth of potentially injured target organisms and which will allow sufficient bacterial growth to allow for sufficient growth of the target organisms in such a way that the pre-enrichment medium will contain at least one viable target organism if such organism existed in the original sample. Examples of suitable non-selective growth media (which are available from Difco Laboratories, Detroit, Michigan) include Tryptic Soy Broth (TSB), Buffered Peptone Water (BPW), Universal Pre-enrichment Broth (UPB), Listeria Enrichment Broth (LEB), and other non-selective or mildly selective media known to those skilled in the art. The preferred pre-enrichment media for use in the present invention is the modified-TSB media.

In the present invention, the pre-enrichment media can also include the addition of a selective inhibitor or mildly selective inhibitor such as an antibiotic. In the method of the present invention, the antibiotic novobiocin can be added to the modified-TSB media prior to the addition and incubation of the sample in order to provide a selection pressure in favor of the growth of *Salmonella*. The novobiocin or sodium salt thereof can be added in a concentration of between 1-50  $\mu\text{g/ml}$ . Preferably, the concentration of the novobiocin is approximately 20-40  $\mu\text{g/ml}$  with the preferred concentration being approximately 25  $\mu\text{g/ml}$  of pre-enrichment media. Alternatively, other antibiotics such as vancomycin, penicillin, ampicillin, and amikacin can be utilized in appropriate concentrations. Other target organism selection enhancing additives known to those skilled in the art can be added to the pre-enrichment broth to favor the growth of the selected

target microorganism over other competing non-target microorganisms present in the sample.

Following the pre-enrichment step or resuscitation of the injured target microorganisms in the sample, a selective enrichment step is carried out to further inhibit competing non-target microorganisms which have not been inhibited during the pre-enrichment step. For example, when *Salmonella* is the target microorganism, the selective enrichment step is carried out to further inhibit other Gram (-) bacteria thereby promoting or encouraging the growth of *Salmonella* in the mixed culture.

At least one selective inhibitor or mixtures thereof can be added at the end of the approximately four hour incubation of the pre-enrichment step. Approximately ten milliliters of the inhibitor mixture is added to the container containing the pre-enrichment media and the sample. This time point is critical since the growth levels of the competing non-target microorganisms will increase to levels that would interfere with obtaining positive reactions to the biochemical analysis unless the growth levels of the competing non-target microorganisms are inhibited. The container with the selective inhibitors is further incubated for an additional approximately six to eight hours (when utilizing semi-automated transfer) or for sixteen to twenty hours (when utilizing manual transfer) at a temperature range of approximately 33°C to 43°C depending on the target microorganism as discussed above. This incubation time window is critical as during this time window, the organisms can be transferred at any time as the ratio of target organisms to non-target microorganisms is suitable for biochemical assaying. Preferably, the container with the selective inhibitor(s) is further incubated for an additional approximately eight hours at approximately 35°C ± 2°C.

The selective inhibitor(s) or mixture thereof can include compounds or reagents which are known to inhibit the growth of non-target microorganisms while encouraging the proliferation or growth of the target microorganisms. For example, if the target microorganism is *Salmonella*, it has been found that a mixture of magnesium chloride, malachite green, and crystal violet can be added to the pre-enrichment media present in the container to yield excellent inhibition of non-target microorganisms. The concentration of the malachite green in the solution ranges from approximately 10mg/l to 50mg/l (0.003%-0.03 w/v), the magnesium chloride concentration ranges from approximately 0.5 g/l to approximately 2 g/l (0.1%-1.0% w/v), and the concentration of the crystal violet ranges from approximately 0.001 g/l to 0.010 g/l (0.0005%-0.001% w/v). The preferred concentration of the inhibitor mixture for *Salmonella* includes 30 mg/l malachite green, 0.005 g/l crystal violet, and approximately 1.0 g/l magnesium chloride. Other inhibitors can include bile salts, sodium deoxycholate, sodium selenite, sodium thiosulfate, sodium tetrathionate, sodium sulphacetamide, mandelic acid, selenite-cysteine tetrathionate, sulphamethazine, brilliant green, malachite green, crystal violet, tergitol 4, sulphadiazine, amikacin, and novobiocin.

After incubation of the sample in the presence of the inhibitor(s) of the non-target microorganisms, an enhanced enrichment step is carried out to eliminate a large proportion, approximately 70-90% or more of the samples which are found to be negative for the presence of the target microorganism. At the end of the approximately twelve hour incubation, (four hours pre-enrichment ± eight hours selective enrichment for semi-automated transfer) or at the end of the approximately twenty hour incubation, (four hours pre-enrichment ± sixteen hours selective enrichment for manual transfer) an aliquot of the sample can be transferred from the original container into separate media receiving vessels wherein specific biochemical assays determinative of the presence of the target microorganism are undertaken. The transfer of the aliquot of the sample can be accomplished either by manually transferring a predetermined amount of the sample into the media receiving vessel or can be accomplished by automated means such as that described in a presently pending application United States Serial Number 08/503,081 to Eden et al., filed July 14, 1995, and assigned to the assignee of the present invention.

Alternatively, the sample can be allowed to incubate overnight (approximately sixteen hours) and then can be either manually or semi-automatically transferred from the original container into separate media receiving vessels wherein specific biochemical assays determinative of the presence of the target microorganism are undertaken. The transfer of the aliquot of the sample can be accomplished either by manually transferring a predetermined amount of the sample into the media receiving vessel or can be accomplished by automated means such as that described above.

Thus, the present invention provides a system in the form of a manual, semiautomatic or automatic system for rapidly detecting a target microorganism in a sample which may contain both target microorganisms and non-target microorganisms. The system can be generally defined as including a container for containing a sample therein, growth inhibitors for inhibiting the growth of non-target microorganisms, biochemical reagents for identifying the target microorganisms, and detecting assays for identifying the target microorganisms.

The biochemical tests or reagents specific for identification of the target microorganisms can include antibiotics, dyes, and other biochemical reagents indicative of particular microorganisms such as by targeting sugar fermentation, decarboxylation, cleavage by unique enzymes, and/or use of unique combinations of dyes (including fluorescent dyes). Additionally, other reagents used in the detection and identification of microorganisms known to those skilled in the art may be practiced with the present invention. For example, for the detection of *Salmonella*, the aliquot of the sample from the container can be transferred into a media receiving vessel, such as a test tube, which can contain reagents necessary for performing a MUCAP (methylumbelliferyl caprylate) test. Additional tubes may be utilized which contain reagents for further biochemical tests such as H<sub>2</sub>S tests, media specifically designed to indicate the presence of lysine

decarboxylase, ornithine decarboxylase, or arginine decarboxylase in a suitable base media such as LICNR (lysine-iron-cystine-neutral red broth, Difco Manual) base media, fermentation reaction media such as for dulcitol, propylene glycol (PG), glucuronic acid (GA), in a peptone base, and a citrate utilization can also be assayed using Simmon's citrate agar (Difco Manual) or citrate medium (magnesium sulfate, ammonium dihydrogen phosphate, sodium citrate, yeast extract, sodium chloride) placed in the media receiving vessels.

The preferred biochemical reagents or tests for *Salmonella* include a tube which includes an H<sub>2</sub>S strip in the MUCAP test and one tube with reagents indicating lysine decarboxylation.

Following the addition of the aliquot of the post-incubation sample into the separate media receiving vessels or tubes containing the biochemical reagents, the samples are incubated for approximately twelve-fourteen hours with twelve hours being the preferred incubation time (a complete enrichment of twenty-four hours).

The results of the biochemical assays, specific for the target microorganisms, allow for the immediate elimination of those samples which are found to be negative for the presence of the target microorganisms. That is, since the biochemical assays or tests have been specifically designated to indicate the presence of the target organism in a mixed culture including target organism and non-target organisms, a negative reaction to these specific biochemical assays indicates the absence of the target microorganisms in the sample and, therefore, allows for the immediate elimination of the sample from consideration.

The detection of a positive or negative test can be done visually, by fluorescence, gas production, and/or by any other suitable detection means known to those skilled in the art.

Following a positive outcome or result after the specific biochemical assays, further microbiological analyses such as immunoassay, DNA probing, serological analysis, plating on selective media, or the like can be undertaken in order to positively identify those samples which are positive for the target microorganism.

By utilizing the method of the present invention, a sample tested which may have contained an initial inoculum of 1 cell/225 ml of target microorganisms, can be enriched to 10<sup>9</sup>-10<sup>1</sup> cfu/ml after pre-enrichment, 10<sup>3</sup>-10<sup>4</sup> cfu/ml following selective enrichment, and greater than 10<sup>6</sup> cfu/ml following enhanced enrichment. At these final concentrations, biochemical assays performed in mixed cultures can be adequately performed to test for the presence of the target microorganisms.

The present invention also includes a kit containing the pre-enrichment media, selective inhibitors, and specific biochemical assay reagents necessary for performing the rapid screening method for the detection and identification of target microorganisms, all of which are described above.

## EXAMPLES

### Heat Injury Methodology:

**Materials:** Tryptic Soy Broth (Difco Lot #800543)  
Tryptic Soy Agar (Difco Lot #66932JA)  
NaCl (Sigma Chemicals 109F/0498)  
Bacto-Peptone (Difco Lot #37680JA)  
TSA+2% NaCl  
TSA  
0.1% Peptone dilution blanks (9.0 ml per tube)  
Water bath (Blue M)

**Procedure:** The water bath (with an agitator) was set at 55°C and was filled to almost full and turned on to warm up for thirty to forty-five minutes with its lid on. Sterile 100 ml capped media bottles each containing 18 ml of TSB were placed into the water bath in a rack for approximately fifteen-twenty minutes prior to the start of the heating assay.

### Enumeration of Cultures Prior to Heating Assay (T<sub>0</sub>)

- 1.) *S. enteritidis* and *S. typhimurium* were grown in TSB for eighteen to twenty-four hours at 35 ± 2°C.
- 2.) Serial dilutions to 10<sup>-6</sup> were performed the day of the assay in sterile peptone dilutions blanks.
- 3.) A 100µl suspension from 10<sup>-6</sup> dilution was spread plated onto the TSA plates. The plates were dried in a biohood for fifteen to twenty minutes and incubated at 35 ± 2°C overnight.

Heat Injury Assay

1.) 2 ml each of *S. enteritidis* and *S. typhimurium* were added to the pre-heated TSB broth in two separate media bottles (1:10 dilution) and were placed in the heated water bath for fifteen minutes.

2.) After fifteen minutes, the media bottles with the inoculum were quickly cooled in an ice bath to stop further injury.

3.) 1 ml aliquots were removed from the media bottles and transferred aseptically to a 9 ml peptone/saline dilution blank ( $10^{-2}$ ).

4.) Serial dilutions up to  $10^{-5}$  were performed.

5.) A 100  $\mu$ l suspension from -2 and -3 dilutions were spread plated onto the TSA and TSA plus 2% NaCl plates. The plates were dried for fifteen-twenty minutes in the biohood and then incubated at  $35 \pm 2^\circ\text{C}$  overnight.

The formulas used for calculating heat injury were as follows:

$$\% \text{ Heat Injury} = 100 \times \frac{\text{cfu/ml on TSA after 15 minutes} - \text{cfu/ml on TSA} + 2\% \text{NaCl after 15 minutes}}{\text{cfu/ml on TSA after 15 minutes}}$$

$$\% \text{ Heat Kill} = \frac{\text{cfu/ml on TSA @ } T_0 - \text{cfu/ml on TSA @ } T_{15}}{\text{cfu/ml on TSA @ } T_0} \times 100$$

Results:

Utilizing the method set forth above for heat injuring bacterial cells to simulate target cells which are sub-lethally injured during testing of a sample, it was found that the injury times in a  $55^\circ\text{C}$  water bath necessary to create a 92-98% injury population ranged from fifteen to thirty minutes.

Growth Medium Selection:

In order to ascertain a suitable growth medium for the resuscitation and growth of microorganisms in a sample, modified-TSB and REVIVE media were compared for their performance with *S. typhimurium*. The six hour plate count before transfer of the inocula to a Biotek light absorbance reader (Bio-Tek Instruments, Winooski, Vermont) is shown in Table 3. A 4. 1 cfu/ml inocula of heat injured *S. typhimurium* was inoculated into modified-TSB, REVIVE, and lactose broth at approximately  $35^\circ\text{C} \pm 2^\circ\text{C}$ . The results of the growth of the injured *S. typhimurium* in the various growth media is shown in Figure 1. Modified-TSB was found to be comparable to the REVIVE media using the Biotek apparatus. Referring to Table 3, the sample grown in the modified-TSB media was found to grow faster and have more cell mass than the sample grown in the REVIVE medium after six hours of incubation indicating greater resuscitation.

Growth Inhibitors of Non-Target Microorganisms:

Several inhibitors were tested with three *Salmonella* (2 *enteritidis*, 1 *typhimurium*) and six competing organisms ((4 Gram negative (-)), 2 Gram positive (+)). All organisms were tested both with and without the inhibitor, individually, and in both modified-TSB and lactose broth (LB). Growth controls were tested without inhibitors for modified-TSB and LB as shown in Figures 2 and 3, respectively.

The organisms tested were *S. enteritidis muenster* (Sm), *S. enteritidis heidelberg* (Sh), *S. typhimurium* (St) and competing microorganisms *Enterobacter aerogenes* (Ea), *E. coli* (Ec), *Citrobacter freundii* (Cf), *Proteus vulgaris* (Pv), *E. faecalis* (Ef), and *S. aureus* (Sa).

The first inhibitor tested was deoxycholate/Brilliant Green 0.25%/0.001%, respectively, the results shown in Figures 4 and 5. Figure 3 illustrates the results of modified-TSB having no deoxycholate/Brilliant Green added thereto. The *Salmonella* species and the Ea grew well in the modified-TSB media with the deoxycholate/Brilliant Green inhibitor. The Ec growth was delayed. The growth of Cf, Pv, Ef, and Sa, were inhibited as shown in Figure 4. Referring to Figure 5, the growth of the *Salmonella* species and Ea in the LB with deoxycholate/Brilliant Green added thereto is illustrated in Figure 5. The growth of Ec was delayed and the growth of Cf, Pv, Ef, and Sa, were inhibited.



Referring to Figure 6, the growth of the test organisms set forth above is illustrated for modified-TSB having the inhibitor combination novobiocin/Brilliant Green (0.0001%/0.001%), respectively, added thereto. The growth of Sm was delayed and the growth of all other test organisms was inhibited. The identical experiment was performed using LB with novobiocin/Brilliant Green and the growth of all test organisms was inhibited.

The effect of the inhibitor tergitol 4 (XLT4 Supplement, Difco Laboratories, Detroit, Michigan) at a concentration of approximately 0.046% was tested in both modified-TSB and LB in order to determine its effect on the test organisms. As shown in Figure 7, in modified-TSB the tergitol 4 inhibited Sa growth. The growth of Pv was delayed and the other test organisms grew well. Referring to Figure 8, the effect of the tergitol 4 in LB illustrates that the growth of Sa and Pv were inhibited while the other test organisms grew well.

The effect of the inhibitor combination tergitol 4/Brilliant Green, 0.046%/0.001%, respectively, in modified-TSB is illustrated in Figure 9. The *Salmonella* species, Ea, and Ec, grew well in the presence of this inhibitor combination. The growth of the Cf was slightly delayed, and the growth of Pv, Ef, and Sa were inhibited.

The inhibitor combination tergitol 4/Brilliant Green in LB inhibited the growth of all of the test organisms.

As shown in Figures 10 and 11 and Tables 4 and 5, the results of experiments designed to determine both the ideal concentration of and time of administration of novobiocin in modified-TSB are illustrated. The effects of different concentrations of novobiocin ranging from 3.125-50µg/ml on heat injured *S. typhimurium* and *S. brandenberg* were studied. Additionally, studies were also performed on combinations of *Salmonella* and non-*Salmonella* (1-10 cfu/ml of *Salmonella* vs. 10<sup>4</sup> cfu/ml non-*Salmonella*) in order to determine the optimum concentration of novobiocin. It was found that at a concentration of 50µg/ml that most gram positive organisms were inhibited and that most gram negative organisms with the exception of *C. freundii* and *E. coli* were inhibited. However, at this concentration, *Salmonella* were also inhibited. A concentration of approximately 25µg/ml in modified-TSB media was found to be the concentration which inhibited gram positives and some gram negatives, while not inhibiting the growth and resuscitation of *Salmonella*. It was also found that the time of addition of the novobiocin was best performed at the very beginning of the test.

#### Biochemical Reagents-Selective Enrichment:

Biochemical reaction formulations were tested using nineteen *Salmonella* serotypes and fourteen non-*Salmonella* competitors as shown in Table 6. A total of eight reactions were tested with the organisms shown in Table 6. These reactions can be divided into four separate groups:

##### 1.) H<sub>2</sub>S Reaction

- dry strip (H<sub>2</sub>S-S) vs. media (H<sub>2</sub>S-M)

##### 2.) Fermentation Reactions

- dulcitol (DUL)
- glucuronate (GA)
- propylene glycol (PG)

##### 3.) Decarboxylase Reaction

- lysine (LYS)
  - ornithine (ORN)
- } in modified LICNR media

- arginine

##### 4.) Enzyme Cleavage Reactions

###### A. Chromogenic (Bromochloroindolyl) -BCI

- glucuronide
- caprylate

###### B. Fluorogenic (Methylumbelliferyl)-MU

- glucuronide (MUG)
- caprylate (MUCAP)

A two tube system can be utilized wherein a first tube (#1) can contain:

H<sub>2</sub>S- strip (black) or  
caprylate (fluorescence)

and a second tube (#2) can contain:

dulcitol (yellow) or  
glucuronide (fluorescence) or  
dulcitol (yellow).

These reactions were investigated as targets for presumptive identification of *Salmonella*. These reactions were chosen because they mostly yield positive test results for *Salmonella* species. Referring to Table 6, the numeric number is the sum of positive test results for each organism for the test. Most *Salmonella* tested positive for all tests. The highest scoring competitors were the two *Citrobacters*, with five positive results, however, it is noted that lysine is negative for these two samples thereby allowing for a positive distinction between the *Salmonella* species and the *Citrobacter spp.*

#### Example 1

Combined enrichment of *Salmonella Spp.* in competition with Gram negative enteric bacteria. This example illustrates an experiment performed through the enrichment step to the enhanced enrichment step. This test was also performed in order to study the effect of the addition of crystal violet (cv) to 0.1% MgCl<sub>2</sub> + 0.003% malachite green on the reaction media.

#### Materials:

For heat injury of the test organisms, refer to the method and material set forth above.

modified-TSB  
(Difco, Lot #91198)  
modified-TSB media (Difco)  
sodium novobiocin (Difco R<sub>x</sub> P5670-42B)  
magnesium chloride, anhydrous (Sigma, Lot #86F-3524)  
malachite green, oxalate salt (Sigma, Lot #122H0257)  
Bacto-crystal violet (Difco, Lot #727317)  
XLD-Agar (Difco, Lot#85746JD)  
XLD-Agar Plates  
Stomacher bags (Seward Stomacher Bags, Seward Medical)  
Bacterial Strains Used: *S. typhimurium* ATCC 14028, *S. brandenberg* ATCC 6955, *E. aerogenes* ATCC 13048, *E. coli* ATCC 25922, and *H. alvei* ATCC 23280.

#### Methods:

See heat injury methodology set forth above.

#### Pre-Enrichment Assay:

Heat injured *S. typhimurium* and *S. brandenberg* were serially diluted in a 1:10 dilution to a concentration of 10<sup>-2</sup>. Two stomacher bags each containing 225 ml of modified-TSB + novobiocin (25 µg/ml) was inoculated with either the *S. typhimurium* or *S. brandenberg* into two separate bags.

The competitive microflora was prepared and added as follows: a culture of each competitor in 0.5 ml McFarland's (approximately 1 x 10<sup>8</sup> cfu/ml) was serially diluted to achieve a concentration of approximately 10<sup>5</sup>. Approximately 33µl of each competitor organism was added to the stomacher bag containing the *Salmonella* and media to achieve a final concentration of competitors of approximately 1 x 10<sup>4</sup> cfu/225 ml. Six stomacher bags were analyzed and arranged as follows:

## EP 0 877 092 A1

(1) *Salmonella typhimurium* plus competitors with crystal violet, (2) *Salmonella brandenberg* plus competitors with crystal violet, (3) competitors alone with crystal violet, (4) *Salmonella typhimurium* plus competitors without crystal violet, (5) *Salmonella brandenberg* plus competitors without crystal violet, and (6) competitors alone without crystal violet.

All of the bags were mixed thoroughly and incubated at approximately 35+/-2°C for approximately four hours. After four hours of incubation, selective inhibitors were added as follows: into all bags to contain crystal violet, 1 ml of a 100 X stock of 0.1% crystal violet was added along with 10 ml of 20 X stock of 0.0003% malachite green + 0.1 % MgCl<sub>2</sub>. 10 ml of the 20 X stock of 0.003% malachite green + 0.1% MgCl<sub>2</sub> was added to all bags marked without crystal violet. The bags were then incubated for an additional eight hours after which the samples were transferred to biochemical reaction reagents as described below.

After transferring the samples to reaction media tubes, the bags were all incubated for an additional twelve hours. Following this incubation, (twenty-four hours from start) plate counts were performed on TSA and XLD plates.

### Results:

#### Initial Inocula:

*S. typhimurium* = 610 c.f.u./225ml (92% heat injury)  
*S. brandenberg* = 101 c.f.u./225ml (91% heat injury)  
Others (*E. coli*, *E. aerogenes*, and *H. alvei*) = 1.0 X 10<sup>4</sup> c.f.u./225ml.

#### Twenty-two Hour Plate Counts: with crystal violet

##### *S. typhimurium*/competitors

*Salmonella* = 2.2 x 10<sup>8</sup> cfu/ml

Competitors = 5.7 x 10<sup>8</sup> cfu/ml

##### *S. brandenberg*/competitors

*Salmonella* = 3.1 x 10<sup>8</sup> cfu/ml

Competitors = 7.0 x 10<sup>8</sup> cfu/ml

### Example 2

#### Complete enrichment model of Salmonella with *E. coli*, *E. aerogenes*, and *Hafnia alvei*.

The purpose of this experiment was to determine if a dye/inhibitor system used during both the pre-enrichment stage and selective enrichment stage would have an effect on the color and fluorescence reactions utilized as negative *Salmonella* screens. Additionally, this experiment was performed in order to obtain population and growth information from a true model system.

### Materials:

disodium L-cystine (Difco, P-2401-20B)  
ferric ammonium citrate (Sigma, 67F-0825)  
salicin (Difco, #798545)  
sodium thiosulfate (Sigma, S-1648)  
mannitol (Difco, 7497L)  
dextrose (Difco, 709544)  
yeast extract (Difco, 0127-19-9)  
L-lysine HCl (Difco, P-0705-20B)

tryptone, pancreatic digest of casein #3 P-5940)  
 neutral red (Sigma, N-7005)  
 H<sub>2</sub>S strips (Difco, 1626-30-6)  
 MUCAP (Research Organics, 0179 M-2)  
 5 phenol red (Difco R<sub>x</sub>, 5533N)  
 dulcitol (Sigma, D-0256)  
 ethyl alcohol, anhydrous (Spectrum Chemical)  
 sodium hydroxide  
 hydrochloric acid  
 10 caprylic acid (Sigma, C-5038)  
 sodium citrate (Difco, P-7260-20B)  
 magnesium sulfate (Sigma, M-1880)  
 ammonium dihydrogen phosphate (Sigma, A-1645)  
 dipotassium phosphate (Difco, P-6500-20B)  
 15 sodium chloride (Difco, 7240 R<sub>x</sub> 6081R)  
 magnesium sulfate (Sigma, M0250)  
 ferrous sulfate (Sigma, F-7002)

#### Reaction Media:

##### 1.) MUCAP + H<sub>2</sub>S Strip Media-citrate based

Ingredient	g/l
magnesium sulfate	0.2
ammonium dihydrogen phosphate	1
dipotassium phosphate	1
sodium chloride	5
yeast extract	0.3
sodium citrate	5
sodium thiosulfate	0.1
pH of media	7.83

##### 2.) Modified LICNR Media

Ingredient	g/l
yeast extract	3
tryptone	5
L-lysine	10
L-cystine	0.1
mannitol	5
dextrose	1
salicin	1
neutral red	0.025
pH of media	6.2+/-0.2°C @ 25°C

Both media 1 and 2 were prepared and dispensed in 9 ml aliquots in media tubes and autoclaved for fifteen minutes at 121°C.

A total of six bags were prepared and contained the modified-TSB media including Novobiocin as set forth above. The bags were inoculated with bacteria as follows:

Bag 1 - St, Ec, Ea, Ha  
 Bag 2 - Sb, Ec, Ea, Ha  
 Bag 3 - Ec, Ea, Ha

Bag 4 - St, Ec, Ea, Ha

Bag 5 - Sb, Ec, Ea, Ha

Bag 6 - Ec, Ea, Ha

5 Bags 1-3 included crystal violet. Bags 4-6 had no crystal violet added as described in the previous example.  
 After twelve hours of incubation including both pre-enrichment and selective enrichment, 1 ml from each bag was transferred into two reaction tubes containing: (1) MUCAP/H<sub>2</sub>S and (2) Lys decarboxylase reagents. The H<sub>2</sub>S strip was folded over the top of the media tube and sealed down with a cap. Usually, half of the length of the strip was used for each sample. The H<sub>2</sub>S strip was incorporated with the MUCAP tube only. The tubes were then incubated to facilitate the reaction for a further twelve hours. Following the twelve hour incubation, the results of the reactions were observed and a count of the bacteria present in the reaction tubes was performed.

A time line for this example is as follows:

10 A.M. inoculate bags with injured *Salmonella* plus competitors  
 15 2 P.M. add inhibitors  
 10 P.M. transfer to biochemical reactions  
 10 A.M. read biochemical reactions for positive results  
 Note: The reactions were read at 8:30 A.M. = 22.5 hours for this experiment and at twenty-four hours for Lys.

#### 20 Procedure For MUCAP Test Results

- 1 ) The samples were incubated in citrate broth with sodium thiosulfate for ten to twelve hours;
- 2 ) A 1 ml aliquot from the sample was transferred to a separate glass reaction tube;
- 25 3 ) A 1% MUCAP stock solution was prepared;
- 4 ) A 1 µl loopful of stock MUCAP solution was added to the 1 ml aliquot of sample; and
- 30 5 ) The fluorescence reaction was allowed to develop and was read every five minutes to at least fifteen minutes and was observed with a Woods lamp.

The results for this example are set forth in Table 7.

35 The MUCAP reagent can either be added directly to the sample by adding an aliquot of MUCAP solution or can be added by first disposing the MUCAP reagent onto a solid support, such as a paper strip. The concentration of the methylumbelliferyl caprylate in the ethyl acetate can range from approximately 0.05% to 10%. The preferred concentration of the methylumbelliferyl caprylate in the ethyl acetate is approximately 0.1%.

40 A dilution of 50 µl of MUCAP reagent dissolved in ethyl acetate can be added to 5ml of media. Additionally, the reagent can be added in a dry state via an impregnated strip and then the strip can be added to the media in order to introduce the MUCAP reagent as set forth below in Example 3. Results for use in the liquid form is the same as those shown in Table 7. Previously, the MUCAP reagent could not be used in liquid form as it had a tendency to precipitate from solution. However, as used herein, the MUCAP reagent dissolved in ethyl acetate can be used either directly in liquid form or the MUCAP reagent dissolved in ethyl acetate can be dried onto a solid substrate and the utilized, in both cases the precipitation has been eliminated.

TABLE 7

	1	2	3	4	5	6
MUCAP *	+++	+++	-	+++	+++	-
H <sub>2</sub> S	+	+	-	+	+	-
Lys	-	-	-	-	-	-

\* after 15 minutes of development

#### 55 H<sub>2</sub>S Results:

+ = any black precipitate formed on strip  
 - = white strip

Lysine Results:

- + = yellow  
- = remains red

MUCAP Results:

- = media background or equal  
+ = very weak fluorescence above background  
++ = fluorescence stronger, may be determined positive without reference  
+++ = good, bright positive  
++++ = strong fluorescence

The MUCAP results were best observed after at least fifteen minutes of incubation.

H<sub>2</sub>S results:

- St samples were positive after nineteen hours.  
Sb samples were positive after 20.5 hours.

MUCAP + NaOH Results:

After 22.5 hours, 10µl of 0.1N NaOH was added to a 1ml test aliquot of each sample in order to increase the observed fluorescence in the tube. The results of the addition of the sodium hydroxide to the MUCAP media are shown in Table 8.

TABLE 8

MUCAP w/NaOH after fifteen minutes.						
Sample	1	2	3	4	5	6
	++++	++++	+	++++	++++	+

It was found that the fluorescence was increased with very little effect on the negative samples (3 and 6) which contained only competitors.

Lysine Results:

After thirty hours of incubation, the lysine test became differential as shown in Table 9.

TABLE 9

1	2	3	4	5	6
+ yel	+ yel	- red	+ yel	+ yel	- red

Bacterial counts were taken at 22.5 hours in the MUCAP/H<sub>2</sub>S tubes. Samples 1, 2, 4, and 5 were counted on XLD media and samples 3 and 6 were counted on TSA media. The results of these counts is set forth in Table 10.

TABLE 10

	1	2	3	4	5	6
<i>Salmonella</i>	9.5x10 <sup>8</sup>	5.5x10 <sup>8</sup>	-	8.9x10 <sup>8</sup>	9.1x10 <sup>8</sup>	-
Competitors	NT	NT	<10 <sup>5</sup>	NT	NT	7.7x10 <sup>8</sup>
NT: not tested						

It was noted that the competitors in sample 3 (bag 3) were reduced to less than 10<sup>5</sup> cfu/ml from the enrichment procedure.

As shown in Table 2, as compared with conventional methods of *Salmonella* testing, the method of the present invention described above allows for the rapid detection of target microorganisms in a sample which can contain non-

target organisms as well as target organisms by specifically tailoring and/or modifying biochemical assays so that they can be performed on a potentially mixed culture thereby greatly reducing the time required for obtaining useful test results.

### 5 Example 3

10 In order to evaluate the performance of the rapid microorganism detection system of the present invention, samples were tested in different conditions. First, a sample was tested which was spiked with *S. typhimurium* and a sample without being spiked with the *S. typhimurium*. A four hour pre-enrichment @ 37°C followed by eighteen hours of selective inhibition at 35°C ± 2°C was performed. Lysine, MUCAP, and H<sub>2</sub>S tests were performed for up to six hours @ 42°C for *Salmonella* detection.

### Material

15 Chicken - Tyson  
Shrimp - shell-on  
Red dye - McCormick  
Dried egg yolk - Sigma E-0625  
Eggs - large, AA  
20 novobiocin - Difco 48580JA  
TSA (Difco 97897)  
Hektoen Enteric Agar  
Stomacher Bags  
Stomacher Apparatus - Tekmar Lab Blender  
25 Lysine-6 for Lysine reaction

Ethyl Acetate (Sigma 27, 052-0)  
0.5% peptone for H<sub>2</sub>S reaction  
MUCAP strips  
30 MgCl<sub>2</sub> (anhydrous) - Sigma M-8266  
Malachite green - Sigma M-6880  
Crystal violet - Difco 727317  
*S. typhimurium* ATCC 14028  
TSB (Difco 97214)  
35 0.85% NaCl - Difco 94669  
Modified-TSB - Difco 91198  
H<sub>2</sub>S strips - Difco 95855 JB  
Commercial Bleach - Clorox  
SDS - J.T. Baker (L050-5)  
40 sterile containers  
sterile scalpels  
sterile glassware  
pipets, tips, etc.

### 45 Methods:

#### Day-1

(1) Modified-TSB: 3 liters total:  
50 98.4g of modified-TSB powder was dissolved in 3 liters of water which was gently heated. The mixture was autoclaved for fifteen minutes @ 121°C, 15 lbs pressure and allowed to cool overnight to room temperature.

(2) Selective inhibitors: 200 ml total:  
4.0g of MgCl<sub>2</sub> was placed into 200 ml of distilled water and dissolved completely, 0.12g of malachite green  
55 was added to the MgCl<sub>2</sub> mixture 0.02g of crystal violet was added to the mixture. The mixture was heated gently to dissolve the solids completely. The mixture was allowed to cool to room temperature and was filter sterilized using a 0.22 µm filter.

(3) *S. typhimurium*: A sub-culture was prepared in TSB and incubated overnight @ 37°C. A 100 µl sample of a 10<sup>-3</sup> dilution from an injured culture (injured as described above) was added to each bag labeled "spiked."

(4) MUCAP strips: A 0.1% solution of MUCAP was prepared in ethyl acetate. Blank filter paper strips were dipped into the 0.1% stock solution and were completely air dried.

(5) Lysine-6 for lysine reaction:

Magnesium sulfate - 0.1 g  
 Ammonium phosphate monobasic - 0.5g  
 Potassium phosphate dibasic - 0.5g  
 NaCl - 2.5g  
 Yeast Extract - 0.15g  
 Lysine HCl - 5.0g  
 Brown Cresol Purple - 2.5ml of 0.2% stock  
 Lactose - 1.25g  
 Sucrose - 1.25g

The above were dissolved in 1L of distilled water, pH was adjusted to 5.1 and sterilized at 121 °C and 15 lbs pressure for 15 minutes. Samples were dispensed into sterile 5ml tubes.

(6) 0.5% peptone

Into 100ml D1 H<sub>2</sub>O. 0.5g of peptone was dissolved and autoclave at 121 °C and 15 lbs pressure for 15 minutes. The solution was allowed to cool to room temperature and 5ml samples were dispensed into sterile tubes.

(7) Hektoen Enteric Agar: A one liter volume was prepared according to the manufacturers' instructions and stored overnight at room temperature.

Day-2

(1) Stomacher Bags 225 ml of modified-TSB was dispensed into sterile stomacher bags and labeled as follows:

1. Dried egg yolk (spiked)	6. Dried Egg Yolk (unspiked)
2. Eggs (spiked)	7. Eggs (unspiked)
3. Chicken (spiked)	8. Chicken (unspiked)
4. Shrimp (spiked)	9. Shrimp (unspiked)
5. Dye (spiked)	10. Dye (unspiked)
	11. Control (spiked, no meat)

(2) Eggs:

A one liter wash solution was prepared as follows:

One gram of SDS was added to 992 ml of H<sub>2</sub>O and 8 ml of commercial bleach was added.

The eggs used in the experiment were washed with a stiff brush and allowed to drain. The washed eggs were allowed to soak for thirty minutes in the SDS-bleach solution. The yolks were aseptically removed and twenty-five grams of each yolk was placed into the appropriately labeled stomacher bags. The bags were stomached for approximately two minutes using the Tekmar stomacher. A one milliliter aliquot was removed for pH test purposes. Another one milliliter aliquot was removed for a serial dilution in saline tubes of 10<sup>-2</sup> and 10<sup>-3</sup> dilutions. One milliliter of 10<sup>-2</sup> and 10<sup>-3</sup> were plated into empty sterile petri dishes and approximately twenty milliliters of liquid plate count agar was poured into the Petri dishes and was allowed to solidify. The plates were then incubated for 48 hours @ 35°C ± 2°C.

(3) Chicken:

Two samples of approximately 25g of slashed chicken were aseptically weighed-out. The samples were placed into stomacher bags and were stomached for approximately two minutes. The stomached material was then sampled and plated as described in (2).

(4) Dried Egg Yolk:



Two samples of approximately 25g of dried egg yolk were aseptically weighed-out and placed in stomacher bags. The samples were stomached and aliquots were removed for pH testing purposes and plating purposes as set forth in (2)

(5) Shrimp:

Two samples of approximately 25g of shrimp slashed in half were aseptically weighed-out and transferred into sterile stomacher bags containing modified-TSB. The samples were handled as described in (2) above.

(6) Novobiocin:

Lyophilized novobiocin was removed from storage @ 4°C and allowed to equilibrate to room temperature. Each vial of lyophilized novobiocin contained 20 mg of novobiocin. Sterile distilled water (3.55 ml) was aseptically injected into each vial to obtain a stock solution of 5.63 mg/ml. One milliliter of novobiocin stock solution was aseptically added into each stomacher bag including the stomacher bag labeled as control (bag 11) to achieve a final concentration of 25 µg/ml. Following the addition of the novobiocin into each stomacher bag, the contents of each bag was mixed by vigorous shaking in order to distribute the novobiocin throughout the mixture.

(7) Plate Counts at Time Zero ( $t_0$ ):

100µl samples of inoculum from the  $10^{-2}$  and  $10^{-3}$  dilutions were plated onto Hektoen Enteric Media as well as onto TSA plates. The plates were then dried for fifteen minutes under laminar flow and incubated for forty-eight hours @ 35°C ± 2°C. Table 7 illustrates the results of the plate counts.

(8) Red Dye:

10 ml of dye was added into 225ml modified-TSB. The procedure as outlined in (2) detailed above was then followed.

(9) Pre-enrichment:

The stomacher bags (with samples and control) were incubated @ 35°C ± 2°C for four hours. Plate counts were taken at  $T_4$  and the results are shown in Table 8.

(10) Selective Inhibition:

The eleven bags were removed from incubation @ 35°C ± 2°C after four hours. 10ml of selective inhibitor was added and mixed. The bags were then incubated @ 35°C ± 2°C overnight.

Day-3

(11) Biochemical Reactions:

(a) Lysine:

5ml of Lysine-6 was dispensed into sterile tubes.

50µl (1:100 dilution) from each bag was added to a separate tube containing Lysine-6. All of the tubes, along with a negative control, were incubated @ 42°C for up to six hours. Each tube was read every hour. The results are shown in Table 11.

(b) H<sub>2</sub>S:

5ml of 0.5% peptone was dispensed into sterile tubes.

100µl from each bag was dispensed into separate tubes containing 0.5% peptone.

An H<sub>2</sub>S strip was hung into each tube along with a negative control tube. The tubes were incubated at 42°C for up to six hours. The results are shown in Table 12.

(c) MUCAP:

5ml of Lysine-6 was dispensed into sterile tubes.

100µl of test material from each bag was placed into separate tubes containing the Lysine-6.

A MUCAP strip was hung into each tube including the negative control tube, the tubes were not inverted and were allowed to incubate @ 42°C for up to six hours. The tubes were then momentarily inverted and any changes were read after five, ten, and fifteen minutes at room temperature. The results are shown in Table 13.

(12) Interpretation of Biochemical Reactions:

(a) Lysine: A change from yellow coloration to purple coloration indicated a positive reaction.

(b) H<sub>2</sub>S: Blackening of the H<sub>2</sub>S strip indicated a positive reaction.

(c) MUCAP: Fluorescence at long wavelength U.V. light of 1 + or better indicated a positive reaction.

(13) Plate Counts at T<sub>22</sub>:

(a) T<sub>22</sub> (22 hours post inoculation): Following overnight incubation with the selective inhibitors, material from each stomacher bag was diluted into saline blanks (10<sup>-1</sup>, 10<sup>-2</sup>, and 10<sup>-3</sup>). 50μl of 10<sup>-2</sup> and 10<sup>-3</sup> were plated onto TSA and HE plates. The plates were incubated at 35°C ± 2°C overnight to determine if *S. typhimurium* or any other contaminant had broken through. The results are shown in Table 9 where (TNTC) indicates that the colonies were too numerous to count.

(b) T<sub>27</sub> (27 hours post incubation): After biochemical reaction with lysine had turned positive (if not, after six hours of incubation @ 42°C) a loopful of reaction mixture was plated onto HE plates and incubated overnight at 35°C ± 2°C and was then observed for colony characteristics. The results are shown in Table 10.

RESULTS:

(1) pH of the sample foods in modified-TSB (following stomaching):

Dried Egg Yolk	6.8
Egg Yellows	7.0
Shrimp	7.22
Chicken	6.68
Dye	7.21
Control (modified-TSB only)	7.22

(2) Percent injury to *S. typhimurium* culture:

The percent injury was calculated as previously described and was 90.3%.

(3) The T<sub>0</sub> plate counts from the stomached samples are shown in Table 7.

CONCLUSIONS:

At time zero (T<sub>0</sub>), none of the samples showed any typical colonies on HE agar (see Table 7). However, the chicken and shrimp samples showed 30 c.f.u./ml at 10<sup>-2</sup> and colonies too numerous to count (TNTC) at the 10<sup>-2</sup> dilution on TSA suggesting contamination.

After four hours of pre-enrichment, the following was observed on TSA plates (see Table 8):

Dried Egg Yolk - spiked	20 c.f.u./ml of 10 <sup>-2</sup> dilution
Chicken - spiked	20 c.f.u./ml of 10 <sup>-2</sup> dilution
Shrimp - spiked	220 c.f.u./ml of 10 <sup>-2</sup> dilution
Chicken - unspiked	20 c.f.u./ml of 10 <sup>-2</sup> dilution
Shrimp - unspiked	240 c.f.u./ml of 10 <sup>-2</sup> dilution

None of the other bags exhibited any growth after four hours of incubation.

Referring to Table 9, after a total of twenty-two (22) hours of incubation (T<sub>22</sub>), it was observed that on TSA, colonies too numerous to count were obtained from samples other than unspiked egg yolk and dye. This finding indicated that in both spiked and unspiked samples, organisms had a chance to grow and could break through, i.e., the bacterial concentration reaches a population > 1 × 10<sup>6</sup> cfu/ml. The egg yolks and dye were found to be sterile and remained sterile throughout the testing period. Following plating on HE agar, dried egg yolk-spiked showed blue-green colonies indicating the presence of *Salmonella* (see Table 10). Egg yolk-spiked showed similar results. From chicken-spiked and shrimp-spiked, colonies were isolated that were salmon colored having bile precipitation. These colonies could possibly have been *E. coli*. The dye-spiked sample showed typical *Salmonella* colonies on HE plates. Dried egg yolk-

unspiked showed typical *Salmonella* colonies on HE plates indicating that dried egg yolk probably had some *Salmonella* contamination. Egg yolk-unspiked and dye remained sterile throughout the entire testing period. Chicken-unspiked was found to contain some contaminants, possibly *E. coli* which appeared salmon colored with bile precipitation. Shrimp-unspiked was found to contain many salmon colored colonies along with a few *Salmonella*-type colonies on HE agar. The positive control (spiked) was found to contain many *Salmonella* appearing colonies on HE agar.

Colonies from inocula from biochemical reaction tubes after six hours @ 42°C appeared as following:

Dried Egg Yolk-spiked, egg yolk-spiked, dye-spiked, dried egg yolk-unspiked, and the control samples were found to have blue-green colonies with black centers indicating *Salmonella*. Chicken-spiked and shrimp-spiked showed salmon colored colonies with bile precipitation indicating the presence of contaminant possibly *E. coli*. Egg yolk-unspiked and dye-unspiked showed no growth, indicating sterility. Chicken-unspiked was found to be positive for the presence of a contaminant, possibly *E. coli*. Shrimp-unspiked was found to have a mixture of contaminants as well as *Salmonella*-like colonies. The positive control was found to contain only *Salmonella*-appearing colonies.

#### Biochemical Reactions:

Referring to Table 14, dried egg yolk-spiked was found to be positive for reactions indicative of *Salmonella*. This result was expected and was believed to be a true positive result. Egg yolk-spiked tested positive for the presence of *Salmonella*. Chicken-spiked elicited a positive biochemical profile for the presence of *Salmonella*. Shrimp-spiked illustrated a discrepancy which was due to the fact that both the lysine and MUCAP reactions were negative for the presence of *Salmonella*. Dye-spiked was unable to be determined due to the fact that both the lysine and MUCAP reactions were difficult to read due to the dark red color of the reaction media. The H<sub>2</sub>S test was found to be positive for the presence of *Salmonella*. Dried egg yolk-unspiked tested negative utilizing the biochemical profile, however, on HE agar, *Salmonella*-like colonies were present which could indicate the presence of *Salmonella* or other contaminants. Egg yolk-unspiked was found to produce a true negative biochemical profile. Chicken-unspiked demonstrated a weakly positive biochemical profile. This result coupled with the fact that contaminants were isolated on HE agar indicated that the biochemical profile for chicken-unspiked was a true negative result. Shrimp-unspiked produced a discrepancy due to the fact that the biochemical results were positive indicating contamination with *Salmonella*. Dye-unspiked remained sterile and did not produce any colonies throughout the entire test. The control tubes illustrated a strong positive biochemical profile for the presence of *Salmonella* which was indicative of a true positive result.

These results, including the discrepancies for both shrimp-spiked and shrimp-unspiked, demonstrated that the rapid microorganism detection method of the present invention yielded definitive results within a twenty-eight hour total testing period.

#### EXAMPLE 4

This example illustrates an experiment performed in order to evaluate the performance of the rapid microorganism detection method of the present invention by testing spiked (with *S. typhimurium*) and unspiked chicken, shrimp, non-fat dry milk, whole black pepper, ground pork, ground beef, and milk chocolate.

#### Materials:

non-fat dry milk - Saco Mix'n Drink

whole black pepper - McCormick

ground pork sausage

ground beef

milk chocolate - Hershey's Milk Chocolate Bar

TSI Agar - Difco) 72853 JG

The remaining materials were as listed in Example 3.

#### Methods:

The methods were carried out as described in Example 3.

Day-2

(1) Stomacher Bags: 225 ml of modified-TSB w/o minerals was dispensed into sterile stomacher bags and labeled as follows:

Spiked	Unspiked
1. Chicken	8. Chicken
2. Shrimp	9. Shrimp
3. Dry Milk	10. Dry Milk
4. Whole Black Peppers	11. Whole Black Peppers
5. Ground Pork Sausage	12. Ground Pork Sausage
6. Ground Beef	13. Ground Beef
7. Milk Chocolate	14. Milk Chocolate
	15. Control

(2) 25g of test material was weighted into stomacher bags labeled as in (1).

(3) Novobiocin: As described above in Example 3.

(4) Plate Counts at Time Zero ( $T_0$ ):

100 $\mu$ l of inoculum from the  $10^{-2}$  and  $10^{-3}$  dilutions were plated onto Hektoen Enteric media as well as TSA plates. The plates were then dried for fifteen minutes under laminar flow and incubated for forty-eight hours @  $35^\circ\text{C} \pm 2^\circ\text{C}$ . The results are shown in Table 15.

(5) Pre-enrichment:

The stomacher bags (with samples and control) were incubated @  $35^\circ\text{C} \pm 2^\circ\text{C}$  for four hours.

(6) Selective Inhibition:

The bags were removed from incubation @  $35^\circ\text{C} \pm 2^\circ\text{C}$  after four hours in incubation. 10ml of selective inhibitor was added and mixed. The bags were then placed into the incubator @  $35^\circ\text{C} \pm 2^\circ\text{C}$  for overnight incubation. Plate counts were taken and  $T_4$  and the results are shown in Table 16.

(7) Biochemical Reactions:

(a) Lysine:

Performed as described above. Results are shown in Table 17.

(b)  $\text{H}_2\text{S}$ :

Performed as described above. Results are shown in Table 18.

c) MUCAP:

Performed as described above. Results are shown in Table 20.

(8) Interpretation of Results:

Lysine: Change from yellow color to purple color indicated a positive test result; change from yellow color to grey color indicated a weak positive test result.

$\text{H}_2\text{S}$ : Slight blackening the strips indicated a weak positive; moderate blackening indicated a positive test.

MUCAP: 1 + or higher in fluorescence indicated a positive test result.

(9) Plating on HE and TSI at  $T_{22}$ :

Hektoen Enteric Agar: Blue-green colonies with black centers indicated *Salmonella*; orange colonies with bile precipitation indicated contaminants. The results are shown in Tables 21 and 22.

TSI Slants: Alkaline slant, black butt, and gas production indicated the presence of *Salmonella*; yellow slant and butt, and gas production indicated the presence of *E. coli*; and black slant and butt, some yellow color (acid slant) indicated cross-contamination between *Salmonella* and a contaminant. The results are shown in Tables 21 and 22.

RESULTS:

(1) pH of test material in modified-TSB after two minutes of stomaching:

control	7.3
chicken	6.82
shrimp	7.34
dry milk	6.79
pepper	7.22
pork sausage	6.86
ground beef	6.84
milk chocolate	6.96

(2) Number of organisms per bag (spiked samples only):

100µl of  $10^{-4}$  injured culture: 1 c.f.u. on a TSA plate yielded approximately 1 c.f.u./bag at  $T_0$ .

$T_0$  on HE plates:

50µl of  $10^{-2}$  and  $10^{-3}$  dilution were plated onto HE plates. The results are set forth in Table 15.

CONCLUSIONS:

Referring to the results depicted on Table 19, dry milk, whole black pepper, and milk chocolate (all unspiked), had true negative test results. All of the seven spiked samples had a positive test for *Salmonella*. Chicken, shrimp, ground pork sausage, and ground beef (all unspiked) were all found to test positive for *Salmonella* indicating that these products may have been contaminated with *Salmonella* or *Salmonella*-like organisms that yielded positive biochemical test profiles for the presence of *Salmonella*.

In all of the positive test results for the presence of *Salmonella*, blue-green colonies were observed either at twenty-two hours or twenty-eight hours on Hektoen Enteric Agar, which is indicative of *Salmonella*-like organisms.

TSI Slant results from these colonies confirmed the presence of *Salmonella*. However,  $10^{-2}$  and  $10^{-3}$  dilutions were found to have overcrowding of colonies on HE agar. In some cases the overcrowding masked the presence of *Salmonella*.

These results indicate that even when an inoculum of only one cell or colony forming unit (c.f.u.) of *S. typhimurium* was inoculated into food samples, that definitive test results were obtained in less time than that required by prior art methods.

Throughout this application various publications are referenced by citation and number. Full citations for the publication are listed below. The disclosure of these publications in their entireties are hereby incorporated by reference into this application in order to more fully describe the state of the art to which this invention pertains.

The invention has been described in an illustrative manner, and it is to be understood the terminology used is intended to be in the nature of description rather than of limitation.

Obviously, many modifications and variations of the present invention are possible in light of the above teachings. Therefore, it is to be understood that within the scope of the appended claims, the invention may be practiced otherwise than as specifically described.

TABLE I

SYSTEM NAME	PRE-ENRICHMENT INC. TIME/TEMP	SELECTIVE STEP INC. TIME/TEMP	ADDITIONAL STEP INC. TIME/TEMP	ASSAY FORMAT	TOTAL TIME	SENSITIVITY/ SPECIFICITY
AOAC	LACTOSE BROTH TRYPTONE SOYA BROTH NUTRIENT BROTH 24h/35C	SELENITE CYSTINE BROTH (SC) OR TETRATHIONATE BROTH (TT) 24h/35C	NONE	PLATING MEDIA - BRILLIANT GREEN AGAR HEKTOEN AGAR XLD AGAR 24h/35C BISMUTH SULFITE AGAR 48h/35C	3 DAYS	100%
ISO	BUFFERED PEPTONE WATER (BPW) 18h/35-37C	RAPAPORT- VASSILIADIS BROTH (RV) 8h/42C OR SC BROTH 8h/35C	NONE	PLATING MEDIA - BRILLIANT GREEN AGAR 24h/35C	3 DAYS	100%
SALMONELLA- TEK ORGANON-TEKNIKA CORPORATION	LACTOSE BROTH NUTRIENT BROTH 24 +/- 2h/35C	TT BROTH 18-24h/42C AND SC BROTH 18-24h/35C	POST-ENRICHMENT: TRANSFER FROM SELECTIVE MEDIUM TO M-BROTH WITH 100µg/ml NOVIOBIOTIN 16-20h/35C	ELISA- MONOCLONAL ANTIBODIES APPROX. 2hrs	4 DAYS	96.5% 89.3%
TECRA INTERNATIONAL BIOPRODUCTS, INC.	NON-INHIBITORY BROTH 18-22h/35C	TT BROTH AND SC BROTH 6-8h/35C	POST-ENRICHMENT: TRANSFER FROM SELECTIVE MEDIUM TO M-BROTH 16-20h/35C	ELISA - SALMONELLA VISUAL IMMUNOASSAY APPROX. 2 hrs	3 DAYS	98.6% 96.1%
TECRA UNIQUE INTERNATIONAL BIOPRODUCTS, INC.	RAW FLESH FOODS: MODIFIED BPW WITH TERGITOL-7 OTHER FOODS: MODIFIED BPW 16h/35C	NON- SELECTIVE THROUGH DIPSTICK	NONE	ELISA - DIPSTICK ASSAY APPROX. 6 hrs	2 DAYS	

TABLE 1 (Continued)

SYSTEM NAME	PRE-ENRICHMENT INC. TIME/TEMP	SELECTIVE STEP INC. TIME/TEMP	ADDITIONAL STEP INC. TIME/TEMP	ASSAY FORMAT	TOTAL TIME	SENSITIVITY/ SPECIFICITY
TECRA SALMONELLA IMMUNOCAPTURE INTERNATIONAL BIOPRODUCTS INC.	MODIFIED BPW WITH EXTRA PHOSPHATE BUFFER MIN. 16h/35-37C	IMMUNOENRICH MENT: 1 TUBE MODIFIED BPW 1 TUBE OF M- BROTH INCUBATE WITH DIPSTICKS AT ROOM TEMP FOR A MINIMUM OF 20 MIN	NONE	ELISA- DIPSTICK ASSAY APPROX. 6hrs	2 DAYS	
PATH STIK LUMAC RV	BP WATER OR BP WATER WITH BRILLIANT GREEN 16-24h/37C	RV BROTH 16-24h/42C	POST-ENRICHMENT: BP WATER 6-8h/37C	ELISA- DIPSTICK ASSAY 10 MINUTES	3 DAYS	93.0% 96.4%
1-2 TEST BIOCONTROL SYSTEMS INC.	LOW MICRO LOAD: NON-INHIBITORY BROTH 24h/35C HIGH MICRO LOAD: PRE-ENRICH AS IN BAM 24h/35C	NONE	NONE	IMMUNODIFFUSION- WHITE BAND OF PPT RESULTS FROM IMMOBILIZATION OF MOTILE SALMONELLA THROUGH SEMI-SOLID MEDIUM 16-30h/35C	3 DAYS	85-100% 100%
BIO MERIEUX VIDAS BIO MERIEUX MODIFIED ISO METHOD	BP WATER 18h/35-37C	RV MEDIA 6-8h/42C AND SC BROTH 6-8h/35-37C	POST-ENRICHMENT: TRANSFER FROM SELECTIVE MEDIUM TO M-BROTH 18h/42C	ELFA- ENZYME-LINKED FLUORESCENT IMMUNOASSAY 45 MINUTES	3 DAYS	
MODIFIED BAM METHOD	NON-SELECTIVE MEDIUM 18h/35-37C	TT BROTH 6-8h/42C AND SC BROTH 6-8h/35-37C				

TABLE 1 (Continued)

SYSTEM NAME	PRE-ENRICHMENT INC. TIME/TEMP	SELECTIVE STEP INC. TIME/TEMP	ADDITIONAL STEP INC. TIME/TEMP	ASSAY FORMAT	TOTAL TIME	SENSITIVITY/ SPECIFICITY
GENE TRAK GENE TRAK SYSTEMS	RAW MEATS: LACTOSE BROTH 22-24h/35C OTHER FOODS: BAMIAOAC METHOD 22-24h/35C	SC BROTH 16-18h/35C TT BROTH AND SC BROTH 6h/35C	POST-ENRICHMENT: TRANSFER FROM SELECTIVE MEDIUM TO GRAM-NEGATIVE BROTH 12-18h/35C	NUCLEIC ACID- BASED ASSAY APPROX. 2hrs	RAW MEATS: 4 DAYS  OTHER FOODS: 3 DAYS	98.5% 97.5%
DYNABEAD DYNAL AS	BP WATER 16-20h/37C	IMMUNOMAGNETIC SEPARATION (IMS): SELECTIVE BEADS ARE COMBINED WITH PRE- ENRICHMENT 10 MIN/AMBIENT TEMP	NONE	PLATING MEDIA- XLD AND BGA 24h/37C	3 DAYS	
MICRO-SCREEN NEOGEN CORP METHOD A	REVIVE MEDIUM 2h/37C	SC BROTH 18h/45C	NONE	CHROMATOGRAPHY COMBINED WITH GOLD LABELED IMMUNOSORBENT ASSAY (GLISA) APPROX. < 8h	2 DAYS	
METHOD B	BAM METHOD LACTOSE BROTH 22-26h/37C	SC AND TT BROTH 24h/37C AND RV BROTH 24h/44C		ISOLATION ON SELECTIVE AGAR MEDIA 24h/35C	4 DAYS	



TABLE 2

COMPARISON OF THE CONVENTIONAL SALMONELLA TESTING METHOD TO DIFCO'S RAPID SALMONELLA TEST				
Time	Conventional	Competitor Method - Reveal	Difco Rapid Salmonella (Automated Transfer)	Difco Rapid Salmonella (Manual Transfer)
Day 1	Preenrichment: -inoculate and stomach food sample in non-inhibitory media -incubate 24 hours	Primary enrichment <u>Preenrichment</u> -inoculate and stomach food sample in non-inhibitory media <u>Selective enrichment</u> -add highly selective inhibitors at 4 hours to preenrichment	Primary enrichment <u>Preenrichment</u> -inoculate and stomach food sample in media with mild inhibitor <u>Selective enrichment</u> -add highly selective inhibitors at 4 hours to preenrichment Secondary enrichment: -transfer aliquot to reaction media at 12 hours via automated device	Primary enrichment <u>Preenrichment</u> -inoculate and stomach food sample in media with mild inhibitor <u>Selective enrichment</u> -add highly selective inhibitors at 4 hours to preenrichment Secondary enrichment: -transfer aliquot to reaction media at 24 hours via manual transfer
Day 2	Selective enrichment: -transfer 1ml aliquot of preenrichment to selective broth media -incubate 24 hours	Detection: -immunoassay -all samples need to be tested <b>Move to confirmation of positives</b>	Analyze reactions: -eliminate 70-90% negative samples from detection Detection: -immunoassay <b>Move to confirmation of positives</b>	Analyze reactions: -eliminate 70-90% negative samples from detection Detection: -immunoassay <b>Move to confirmation of positives</b>
Day 3	Detection: -Plate all samples on selective agar from selective enrichment broth to obtain pure cultures -incubate 24 hours			
Day 4	Analyze plates for suspect <i>Salmonella</i> colonies <b>Move to confirmation of positives</b>			

TABLE 3

SIX HOUR PLATE COUNT BEFORE TRANSFER TO BIOTEK		
BROTH	TSA	TSA+NaCl
REVIVE	1.5E+04	1.5E+04
ESP-TSB	2.1E+04	2.1E+04
LACTOSE	8.6E+03	4.7E+03

TABLE 4

GROWTH OVER TIME FOR S. TYPHIMURIUM VARIOUS CONCENTRATIONS OF NOVOBIOCIN						
Time	Control	50ug/ml	25ug/ml	12.5ug/ml	6.25 ug/ml	3.125 ug/ml
10h	1.0E+06	2E+03	2.3E+05	6.0E+05	6.5E+05	1.0E+06
12h	2.4E+08	1.6E+04	2.0E+06	5.4E+06	6.5E+06	1.7E+08
14h	6.9E+08	3E+04	9.1E+07	2.3E+08	4.2E+08	8.8E+08
16h	9.4E+08	1E+05	5.2E+08	7.3E+08	9.2E+08	9.6E+08
18h	1.2E+09	9E+05	8.2E+08	9.4E+08	8.7E+08	9.4E+08

TABLE 5

GROWTH OVER TIME FOR S. BRANDENBURG VARIOUS CONCENTRATIONS OF NOVOBIOCIN						
Time	Control	50 ug/ml	25ug/ml	12.5 ug/ml	6.25 ug/ml	3.125 ug/ml
10h	7.0E+07	1.0E+05	1.8E+06	1.5E+07	5.0E+07	3.7E+07
12h	8.8E+08	3.4E+05	1.4E+08	1.2E+09	7.9E+08	8.9E+08
14h	1.0E+09	1.5E+06	6.8E+00	1.4E+09	1.9E+09	1.9E+09
16h	3.0E+09	3.3E+07	1.3E+09	1.1E+09	1.4E+09	1.2E+09
18h	2.1E+09	7.0E+07	3.2E+09	3.6E+09	5.5E+09	1.0E+09

TABLE 6

## Reaction Results for Salmonella and Competitors

SALMONELLA	DUL	PG	GA	LYS	ORN	H2S-M	H2S-S	CIT	TOTAL
<i>S. agona</i>	+	+	+	+	+	+	+	+	8
<i>S. anatum</i>	+	+	+	+	+	+	+	+	8
<i>S. braenderup</i>	+	+	+	+	+	+	+	+	8
<i>S. brandenburg</i>	+	+	+	+	+	+	+	+	8
<i>S. derby</i> *	+	+	+	+	+	+	+	+	8
<i>S. heidelberg</i>	+	+	+	+	+	+	+	+	8
<i>S. infantis</i>	+	+	+	-	-	+	+	+	6
<i>S. kentucky</i>	+	+	+	+	+	+	+	+	8
<i>S. montevideo</i>	+	+	+	+	-	+	+	+	7
<i>S. muenchen</i>	+	+	+	+	+	+	+	+	8
<i>S. muenster</i>	+	+	+	+	+	+	+	+	8
<i>S. newport</i>	+	+	+	+	+	+	+	+	8
<i>S. oranienberg</i> *	+	+	+	+	+	+	+	+	8
<i>S. pullorum</i>	-	-	+	+	+	+	+	+	6
<i>S. reading</i> *	+	+	+	+	+	+	+	+	8
<i>S. schottmuelleri</i> *	+	+	+	+	+	+	+	+	8
<i>S. st. paul</i> *	+	+	+	+	+	+	+	+	8
<i>S. thompson</i> *	+	+	+	+	+	+	+	+	8
<i>S. typhimurium</i>	+	+	+	+	+	+	+	+	8

\* indicates new organisms tested

COMPETITORS	DUL	PG	GA	LYS	ORN	H2S-M	H2S-S	CIT	TOTAL
<i>Acinebacter calcoaceticus</i>	-	-	-	-	-	-	-	-	0
<i>Alcaligenes fecalis</i>	-	-	-	-	-	-	-	-	0
<i>Citrobacter diversus</i>	+	+	+	-	+	-	+	-	5
<i>Citrobacter freundii</i>	+	+	+	-	-	-	w+	+	5
<i>Enterobacter aerogenes</i>	-	-	+	+	-	-	-	+	3
<i>Enterococcus faecalis</i>	-	-	-	-	-	-	-	-	0
<i>E. coli</i>	-	-	+	-	-	-	-	-	1
<i>Hafnia alvei</i>	-	-	+	-	-	+	-	-	2
<i>Klebsiella pneumoniae</i>	-	-	+	+	-	-	-	+	3
<i>Proteus mirabilis</i>	-	-	-	-	+	+	+	-	3
<i>Proteus vulgaris</i>	-	-	-	-	-	-	-	-	0
<i>Providencia alcalifaciens</i>	-	-	-	-	-	-	-	+	1
<i>Serratia marcescens</i>	-	-	-	-	-	-	-	+	1
<i>Staph aureus</i>	-	-	-	-	-	-	-	-	0

TABLE 7

Sample To	# of Colonies on HE Agar: 100 µl Stomached Samples, 37C	
	10 <sup>-2</sup> ; T-24	10 <sup>-2</sup> ; T-48
DEYolk	0	0
Eyolk	0	0
Chicken	0	0
Shrimp	0	0
Dye	0	0
# of Colonies on PCA: 1 ml Stomached Samples: T48(incubation)		
Sample	1.0 <sup>-2</sup>	10 <sup>-3</sup>
DEYolk	0	0
Eyolk	0	0
Chicken	30	8
Shrimp	TNTC	129
Dye	0	0

Table 8

T-4 After four hours of pre-enrichment									
24 HOUR					48 HOURS				
Media	TSA	TSA	HE		HE	TSA	TSA	HE	HE
Sample	10 <sup>-2</sup>	10 <sup>-3</sup>	10 <sup>-2</sup>		10 <sup>-3</sup>	10 <sup>-2</sup>	10 <sup>-3</sup>	10 <sup>-2</sup>	10 <sup>-3</sup>
DEYolk-S	20 c.f.u./ml	0	0		0	1	0	0	0
EYolk-S	0	0	0		0	0	0	0	0
Chicken-S	20 c.f.u./ml	0	0		0	1	0	0	0
Shrimp-S	220 c.f.u./ml	2	0		0	11	2	0	0
Dye-S	0	0	0		0	0	0	0	0
DEYolk	0	0	0		0	0	0	0	0
EYolk	0	0	0		0	0	0	0	0
Chicken	20 c.f.u./ml	0	0		0	1	0	0	0
Shrimp	12	2	0		0	12	2	0	0
Dye	0	0	0		0	0	0	0	0
Control	0	0	0		0	0	0	0	0

TABLE 9

T-22, 37°C After overnight selective inhibition				
Media	TSA	TSA	HE	HE
Sample	10 <sup>-2</sup>	10 <sup>-3</sup>	10 <sup>-2</sup>	10 <sup>-3</sup>
DEYolk-S	TNTC	TNTC	TNTC S. Typhimurium	TNTC S. Typhimurium

TABLE 9 (continued)

T-22, 37°C After overnight selective inhibition				
Media	TSA	TSA	HE	HE
EYolk-S	TNTC	TNTC	TNTC S. Typhimurium	TNTC S. Typhimurium
Chicken-S	TNTC	TNTC	TNTC <u>E. coli</u> , Few st.	TNTC <u>E. coli</u> , Few st.
Shrimp-S	TNTC	TNTC	TNTC <u>E. coli</u> , Few st	TNTC <u>E. coli</u> , Few st.
Dye-S	TNTC	TNTC	TNTC S. Typhimurium	TNTC S. Typhimurium
DEYolk	TNTC	TNTC	TNTC S. Typhimurium	TNTC S. Typhimurium
E. Yolk	0	0	0	0
Chicken	TNTC	TNTC	TNTC <u>E. coli?</u>	TNTC <u>E. coli?</u>
Shrimp	TNTC	TNTC	TNTC <u>E. coli</u> , Few st.	TNTC <u>E. coli</u> , Few. st.
Dye	0	0	0	0
Control	TNTC	TNTC	TNTC S. Typhimurium	TNTC S. Typhimurium

TABLE 10

Colony Morphology:

5Hrs/42C Lysine Rxn on HE Ager (10 µl/plate)

	Sample	Colony Colour	Possible Organism
SPIKED	DEYolk-S	Blue green	<i>S. typhimurium</i>
	Eyolk-S	Blue green	<i>S. typhimurium</i>
	Chicken-S	Salmon, Bile Precipitation	<i>E. coli</i>
	Shrimp-S	Salmon, Bile Precipitation	<i>E. coli</i>
	Dye-S	Blue Green	<i>S. typhimurium</i>
UNSPIKED	DEYolk	Blue Green	<i>S. typhimurium</i>
	Eyolk	No Growth	---
	Chicken	Salmon, Bile Precipitation	<i>E. coli</i>
	Shrimp	Salmon, Bile Precipitation; Blue Green	<i>E. coli</i> <i>S. typhimurium</i>
	Dye	No Growth	---
	Control	Blue Green	<i>S. typhimurium</i>

TABLE 11

## Lysine Reaction

Media	S1	S2	S3	S4	S5	S6	S7	S8	S9
1Hr	Light Green	Light Green	Yell-Green	Yell-Green	Light Red	Yell-Green	Yell-Green	Yell-Green	Grey Green
2Hr	Light Green	Light Green	Yell-Green	Yell-Green	Light Red	Yell-Green	Yell-Green	Yell-Green	Purple Grey
3Hr	Light Green	Light Green	Yell-Green	Yellow	Light Red	Yell-Green	Yell-Green	Yell-Green	Purple
4Hr	Purple	Purple	Purple	Yellow	Fuchsia	Yell-Green	Yell-Green	Yell-Green	Purple
5Hr	Dark Purple	Dark Purple	Purple	Yellow	Fuchsia	Yell-Green	Yell-Green	Purple Grey	Purple
6Hr	Dark Purple	Dark Purple	Purple	Yellow	Fuchsia	Yell-Green	Yell-Green	Grey Purple	Purple
Over-night	V. Dark Purple	V. Dark Purple	Light Grey	Yellow	Dark Fuchsia	Yell-Green	Yell-Green	Grey	Purple

Media	S10	S11	Negative
1Hr	Light Red	Light Green	Yellow
2Hr	Light Red	Light Green	Yellow
3Hr	Light Red	Green Grey	Yellow
4Hr	Amber	Purple	Yellow
5Hr	Amber	Dark Purple	Yellow
6Hr	Amber	Dark Purple	Yellow
Over-night	Amber	Very Dark Purple	Yellow

TABLE 12

H<sub>2</sub>S Reaction

Media	S1	S2	S3	S4	S5	S6	S7	S8	S9
1 Hr	(-)	(-)	(-)	(-)	(-)	(-)	(-)	(-)	(-)
2 Hr	wk1+	wk1+	(-)	wwk1+	wk1+	(-)	(-)	(-)	wk1+
3 Hr	2+	2+	wk1+	2+	2+	(-)	(-)	wk1+	2+
4 Hr	3+	3+	2+	3+	3+	(-)	(-)	wk1+	3+
5 Hr	4+	4+	3+	3+	4+	(-)	(-)	wk1+	4+
6 Hr	4+	4+	3+	4+	4+	(-)	(-)	1+	4+
Over-night	4+	4+	4+	4+	4+	(-)	(-)	1+	4+

Media	S10	S11	Negative
1Hr	(-)	(-)	(-)
2Hr	(-)	1+	(-)
3Hr	(-)	2+	(-)
4Hr	(-)	3+	(-)
5Hr	(-)	4+	(-)
6Hr	(-)	4+	(-)
Over-night	(-)	4+	(-)



TABLE 13

## MUCAP REACTION

Media	S1	S2	S3	S4	S5	S6	S7	S8	S9
5'	BGr	BGr	BGr-1+	BGr	BGr	BGr	BGr	BGr	1+
10'	1+	2+	2+	BGr	BGr	BGr	BGr	wk1+	2+
15'	2+	2+	3+	BGr	1+	BGr	BGr	1+	2+
4 Hr									
5 Hr									
6 Hr									
Over-night									

Media	S10	S11	Negative
5'	BGr	2+	BGr
10'	BGr	3+	BGr
15'	BGr	4+	BGr ..
4 Hr			
5 Hr			
6 Hr			
Over-night			

TABLE 14

Growth History: 1/15/97 &gt; 1/17/97

	Sample	T0-PCA	T4-TSA	T4-HE	T22-TSA	T22-HE	Aft. Rxn 27 Hrs	Lys	H <sub>2</sub> S	Mcap 10'
SPIKED	DEYolk-S	0	20cfu/ml	0	TNTC	T-St	St	+	+	1+
	EYolk-S	0	0	0	TNTC	T-St	St	+	+	2+
	Chicken-S	8	20cfu/ml	0	TNTC	TNTC Ec, Few St	Ec	+	+	2+
	Shrimp-S	129	220cfu/ml	0	TNTC	TNTC Ec, Few St	Ec	Neg	+	Bgr
	Dye-S	0	0	0	TNTC	T-St	St	Neg*	+	Bgr
UNSPIKED	DEYolk	-	0	0	TNTC	T-St	St	Neg	Neg	Bgr
	EYolk	-	0	0	0	0	0	Neg	Neg	Bgr
	Chicken	-	20cfu/ml	0	TNTC	T-Ec	Ec	Wk+	Wk+	Wk1+
	Shrimp	-	240cfu/ml	0	TNTC	T-Ec, Few St	Ec+St	+	+	2+
	Dye	-	0	0	0	0	0	Neg*	Neg	Bgr
	Control	-	0	0	TNTC	T-St	St	+	+	3+

\* Difficult to interpret colour change due to red coloration of dye

TABLE 15

# of Colonies on both HE Agar and TSA: 50µl Stomached Samples, 37°C				
Sample	HE 10 <sup>-2</sup>	TSA 10 <sup>-2</sup>	HE 10 <sup>-3</sup>	TSA 10 <sup>-3</sup>
Chicken	0	0	0	0
Shrimp	2BG	48	0	5
Dry Milk	0	0	0	0
Whole Black Pepper	0	62	0	25
Ground Pork Sausage	30g 13BG	123	2 or	53
Ground beef	0	0	0	0
Milk Chocolate	0	5	0	0

TABLE 16

# of colonies on HE; TSA at T-4				
Media	TSA	TSA	HE	HE
Sample	10*-2	10*-3	10*-2	10*-3
Chicken - spiked	9	6	2	0
Shrimp - spiked	64	12	23BG	2BG
Dry Milk - spiked	0	0	0	0
Whole Pepper - spiked	TNTC	27	0	0
Ground Pork Sausage - spiked	139	37	65 or 20BG	5 or 5BG
Ground Beef- spiked	6	0	40 or	0
Milk Chocolate spiked	3	0	0	0
Chicken - unspiked	5	0	1 or 3 BG	0
Shrimp - unspiked	130	21	18 BG	2BG
Dry Milk - unspiked	1	0	0	0
Whole Pepper- unspiked	TNTC	24	0	0
Ground Pork Sausage - unspiked	129	50	60 or 22BG	4 or
Ground Beef - unspiked	10	0	2 or 3BG	0
Milk Chocolate - unspiked	5	1	0	0
Control	0	0	0	0
or = Orange BG = Blue Green				

TABLE 17

## Lysine Decarboxylation Test Results

Media	S1	S2	S3	S4	S5	S6	S7	S8	S9
1Hr	Yellow	Yellow	Yellow	Yellow	Yellow	Yellow	Yellow	Yellow	Yellow
2Hr	Yellow	Yellow	Yellow	Grey	Grey	Grey	Grey	Yellow	Yellow
3Hr	Yellow	Yellow	Yellow	Purple Grey	Purple Grey	Light Purple	Purple Grey	Yellow	Grey
4Hr	Light Grey	Light Grey	Yellow	Light Purple	Dark Purple	Purple	Purple	Grey	Dark Grey
5Hr	Grey	Grey	Yellow	Dark Purple	Dark Purple	Purple Grey	Dark Purple	Purple Grey	Grey
6Hr	Purple Grey	Grey	Yellow	Dark Purple	Dark Purple	Purple Grey	Dark Purple	Purple	Grey
Over- night	Light Grey	Yellow Grey	Yellow	Dark Purple	Light Grey- Yellow	Light Grey	Dark Purple	Yellow Grey	Yellow Grey

Media	S10	S11	S12	S13	S14	S15	Negative
1 Hr	Yellow	Yellow	Yellow	Yellow	Yellow	Yellow Green	Yellow
2 Hr	Yellow	Yellow	Grey	Yellow	Yellow	Grey	Yellow
3 Hr	Yellow	Yellow	Light Purple	Purple Grey	Yellow	Purple Grey	Yellow
4 Hr	Yellow	Yellow	Dark Purple	Purple	Yellow	Purple	Yellow
5 Hr	Yellow	Yellow	Dark Purple	Purple Grey	Yellow	Dark Purple	Yellow
6 Hr	Yellow	Yellow	Dark Purple	Purple Grey	Yellow	Dark Purple	Yellow
Overnight	Yellow	Yellow	Yellow Grey	Light Purple Grey	yellow	Dark Purple	Yellow

TABLE 18

H<sub>2</sub>S Test Results

Media	S1	S2	S3	S4	S5	S6	S7	S8	S9
1 Hr	Neg.	Neg.	Neg.	Neg.	Neg.	Neg.	Neg.	Neg.	Neg.
2 Hr	Neg.	Neg.	wk +	wk +	neg.	wk +	1 +	wk +	wk +
3 Hr	Neg.	Neg.	1 +	1 +	wk +	wk2+	wk2+	wk+	wk1+
4 Hr	wk +	wk1+	1 +	2+	2+	2+	3+	1+	2+
5 Hr	wk+	1+	wk2+	3+	3+	3+	4+	1+	3+
6 Hr	wk+	2-3+	2+	4+	4+	4+	4+	2+	4+
Over-night	1+	3+	4+	4+	4+	4+	4+	3+	4+

Media	S10	S11	S12	S13	S14	S15	Negative	S15-Retest
1Hr	Neg.	Neg.	Neg.	Neg.	Neg.	Neg.	Neg.	Neg.
2 Hr	+	Neg.	Neg.	wk+	Neg.	Neg.	Neg.	wk1+
3 Hr	1+	Neg.	wk+	2+	Neg.	Neg.	Neg.	1+
4 Hr	1+	Neg.	1+	3+	Neg.	Neg.	Neg.	2+
5 Hr	1+	Neg.	2+	4+	Neg.	Neg.	Neg.	not read
6 Hr	2+	Neg.	4 +	4 +	Neg.	Neg.	Neg.	not read
Over-night	2+	Neg.	4+	4+	Neg.	Neg.	Neg.	4+

TABLE 19

Sample	T0-TSA From 10 <sup>-2</sup>	T0 HE BG Col. From 10 <sup>-2</sup>	T0 PCA From 10 <sup>-2</sup>	T4 HE BG Col.	T22-HE	Tai:T22Hrs	Aft. Ron 29-30 Hrs: HE	TSl:T29- 30 Hrs	Lys Ron	H2S Rxn	Muca p Rxn	Result
Chicken-S	0	0	4.2X10 <sup>3</sup> No Flr.	1X10 <sup>3</sup>	TNTC Or+BG	A+G: Ec	TNTC Or+Bg	A+G: Ec	Wk+	Wk+	Pos	+St
Shrimp-S	4.8X10 <sup>3</sup>	2X10 <sup>2</sup>	TNTC, No Flr.	2.3X10 <sup>3</sup>	TNTC Or+BG	H2S+G: Sl	TNTC Or	A+G: Ec	Wk+	Wk+	Pos	+St
DryMilk-S	0	0	3, No Flr. 3X10 <sup>4</sup>	0	TNTC Or+BG	A+G: Ec	TNTC Or+BG	H2S+G: Sl	Neg	Pos	Pos	+St
WhiBlkPprs-S	6.2X10 <sup>3</sup>	0	TNTC, Flr.	0	TNTC BG	H2S+G: Sl	TNTC BG	H2S+G	Pos	Pos	Pos	+St
GrPctSau-S	1.2X10 <sup>4</sup>	1.3X10 <sup>3</sup>	TNTC, Flr.	2X10 <sup>3</sup>	TNTC Or+BG	A+G Ec	TNTC Or+BG	112S+G: Sl	Pos	Pos	Pos	+St
GrBeef-S	0	0	5.5X10 <sup>3</sup> W K Flr.	0	TNTC Or+BG	H2S+G: Sl	TNTC Or+BG	H2S+G: Sl	Pos	Pos	Pos	+St
MilkChoc-S	0	0	4.8X10 <sup>3</sup> Flr.	0	TNTC Or+BG	H2S+G: Sl	TNTC BG	H2S+G: Sl	Pos	Pos	Pos	+St
Chicken	NA	NA	NA	3X10 <sup>2</sup>	TNTC Or+BG	H2S+G: Sl	TNTC Or+BG	A+G: Ec	Pos	Pos	Pos	+St
Shrimp	NA	NA	NA	1.8X10 <sup>3</sup>	TNTC Or+BG	H2S+G: Sl	TNTC Or+BG	A+G: Ec	Wk+	Pos	Pos	+St
DryMilk	NA	NA	NA	0	TNTC Or	A+G: Ec	TNTC Or	A+G: Ec	Neg	Wk+	Bgr	Neg
WhiBlkPprs	NA	NA	NA	0	0	0	0	0	Neg	Neg	Bgr	Neg
GrPctSau	NA	NA	NA	2.2X10 <sup>3</sup>	TNTC Or+BG	H2S+G: Sl	TNTC Or+BG	A+G: Ec	Pos	Pos	Pos	+St
GrBeef	NA	NA	NA	3X10 <sup>2</sup>	TNTC Or+BG	H2S+G: Sl	TNTC Or+BG	H2S+G: Sl	Pos	Pos	Pos	+St
MilkChoc	NA	NA	NA	0	0	0	0	0	Neg	Neg	Bgr	Neg
Control-S	NA	NA	NA	0	TNTC BG	H2S+G: Sl	TNTC BG	H2S+G: Sl	Pos	Pos	Pos	+St

Flr = Fluorescence; BG = Blue Green Colonies; Or = Orange Colonies; TNTC = Too Numerous to Count; NA = Not Available  
H2S = H2S production & blackening; G = Gas production; A = Acid slant (yellow); Bgr = Background fluorescence

TABLE 20

## MUCAP Test Results

Media	S1	S2	S3	S4	S5	S6	S7	S8	S9
	4+	3+	2+	4+	4+	4+	4+	4+	3+

Media	S10	S11	S12	S13	S14	S15	Neg.
	BGr	BGr	4+	4+	BGr	4+	Neg BGr

TABLE 21

# of colonies on HE and TSI @ T <sub>22</sub>						
Media	TSA	TSA	HE	HE	TSI	TSI
Sample	10 <sup>-2</sup>	10 <sup>-3</sup>	10 <sup>-2</sup>	10 <sup>-3</sup>	10 <sup>-2</sup>	10 <sup>-3</sup>
Chicken-spiked	TNTC	Or+BG	TNTC	TNTC	A + G	EC?
Shrimp-spiked	TNTC	Or+BG	TNTC	TNTC	H <sub>2</sub> S + G	St?
Dry Milk-spiked	TNTC	Or+BG	TNTC	TNTC	A + G	EC?
Whole Pepper-spiked	TNTC	+ BG	TNTC	TNTC	H <sub>2</sub> S + G	St?
Ground Pork Sausage-spiked	TNTC	Or + BG	TNTC	TNTC	A + G	EC?
Ground Beef-spiked	TNTC	Or + BG	TNTC	TNTC	H <sub>2</sub> S + G	St?
MilkChocolate-Spiked	TNTC	Or + BG	TNTC	TNTC	H <sub>2</sub> S + G	St?
Chicken-unspiked	TNTC	Or + BG	TNTC	TNTC	H <sub>2</sub> S + G	St?
Shrimp-unspiked	TNTC	Or + BG	TNTC	TNTC	H <sub>2</sub> S + G	St?
Dry Milk-unspiked	TNTC	Or	TNTC	TNTC	A + G	EC?
WholePepper-unspiked	0	0	0	0	0	0
Ground Pork Sausage-unspiked	TNTC	Or + BG	TNTC	TNTC	H <sub>2</sub> S + G	St?
Ground Beef-unspiked	TNTC	Or + BG	TNTC	TNTC	H <sub>2</sub> S + G	St?
Milk Chocolate-unspiked	0	0	0	0	0	0
Control-S	TNTC	BG	TNTC	TNTC	H <sub>2</sub> S + G	St?
A = acid slant      G = gas      H <sub>2</sub> S = H <sub>2</sub> S prod. blackening EC = <i>E. coli</i> St - <i>S. typhimurium</i>						

TABLE 22

Colony Morphology:			
After 5hrs/42°C Biochemical Rxn on HE Agar (10μl/plate)			
Sample	HE Morphology	TSI Result	Possible Organism
Chicken-spiked	Orange-TNTC/BP BG - few	A + G	<i>E. coli</i> (Hard to pick colonies)

TABLE 22 (continued)

Colony Morphology:			
After 5hrs/42°C Biochemical Rxn on HE Agar (10μl/plate)			
Sample	HE Morphology	TSI Result	Possible Organism
Shrimp-spiked	Orange-TNTC/BP	A + G	<i>E. coli</i> (Hard to pick colonies)
Dry Milk-spiked	Orange-TNTC BP BG - Few	H <sub>2</sub> S + G	St
Whole Black Pepper - spiked	BG-TNTC	H <sub>2</sub> S + G	St
Ground Pork Sausage-spiked	Orange-TNTC/BP BG - Few	H <sub>2</sub> S + G	St
Ground Beef - spiked	Orange-TNTC/BP BG - Few	H <sub>2</sub> S + G	St
Milk Chocolate -spiked	BG-TNTC	H <sub>2</sub> S + G	St
Control-S	BG - TNTC	H <sub>2</sub> S + G	St
Chicken - unspiked	Orange - TNTC/BP BG - Few	A + G	<i>E. Coli</i> (Hard to pick colonies)
Shrimp - unspiked	Orange-TNTC/BP BG - Few	A + G	<i>E. coli</i> (Hard to pick colonies)
Dry Milk - unspiked	Orange-TNTC/BP	A + G	<i>E. coli</i> (Hard to pick colonies)
Whole Black Pepper - unspiked	0	-	-
Ground Pork Sausage - unspiked	Orange-TNTC/BP BG - Few	A + G	<i>E. coli</i> (Hard to pick colonies)
Ground Beef - - unspiked	Orange-TNTC/BP BG - Few	H <sub>2</sub> S + G	St
Milk Chocolate - unspiked	0	-	-
Control	-	-	-
BP = bile precipitation      A = acid slant      G = gas production H <sub>2</sub> S = H <sub>2</sub> S production black color			

#### Claims

1. A rapid screening method for the detection and identification of target microorganisms in a sample which may contain both target microorganisms and competing non-target microorganisms, said method comprising the steps of:
  - (a) pre-enriching a sample in a growth medium;
  - (b) adding at least one inhibitor of said non-target microorganisms to said growth medium to discourage the growth of said non-target microorganisms and to encourage the growth of said target microorganisms;
  - (c) incubating said sample in said growth medium including said at least one inhibitor for a predetermined amount of time;
  - (d) performing biochemical assays specific for identification of said target microorganisms; and
  - (e) detecting the presence of said target microorganisms in said sample.
2. A method according to claim 1, wherein said growth medium is mildly inhibitory.
3. A method according to claims 1 or 2, which further includes the step of screening out samples found to be negative for the target microorganism.



4. A method according to any of claims 1 to 3, wherein the sample includes a mixed culture of microorganisms.
5. A method according to any of claims 1 to 4, wherein step (a) is approximately two to six hours in duration, preferably about four hours in duration.
6. A method according to any of claims 1 to 5, wherein said predetermined time in step (c) is approximately six to eight hours, preferably approximately eight hours.
7. A method according to any of claims 1 to 6, wherein said sample is incubated in step (c) with said growth medium containing biochemical reagents specific for identification of the target microorganism for approximately twelve to fourteen hours, preferably approximately twelve hours.
8. A method according to any of claims 1 to 7, wherein the sample is incubated in step (c) at a temperature ranging from approximately 33°C to 43°C, preferably approximately 35°C.
9. A method according to any of claims 1 to 8, which further includes the step of transferring a predetermined amount of the growth medium containing the sample into a least one media receiving vessel containing biochemical reagents specific for the identification of the target microorganisms (which transferring step is typically substantially automated, or substantially manual).
10. A method according to claim 9, wherein said transferring step further includes the step of performing an immunological assay on the sample contained in the predetermined amount of growth medium.
11. A method according to any of claims 1 to 10, wherein said sample is at least partially homogenized in the growth medium.
12. A system for rapid screening for the detection and identification of target microorganisms in a sample which may contain both target microorganisms and competing non-target microorganisms, said system comprising:
  - (a) container means for containing said sample and a growth medium;
  - (b) growth inhibition means for inhibiting the growth of non-target microorganisms in said growth medium to discourage the growth of non-target microorganisms and to encourage the growth of said target microorganisms;
  - (c) biochemical reagent means for identification of said target microorganisms in the sample; and
  - (d) detection means for detecting the presence of said target microorganisms in said sample.
13. A system according to claim 12, or a method according to any of claims 1 to 8, wherein said growth medium includes modified-TSB and an antibiotic, which preferably includes novobiocin.
14. A system according to claim 12 or 13, or a method according to any of claims 1 to 8, wherein said growth inhibition means include a mixture of magnesium chloride, malachite green, and Crystal violet, or includes selenite-cysteine tetrathionate.
15. A system according to any of claims 12 to 14, or a method according to any of claims 1 to 8, wherein said biochemical reagent means include an H<sub>2</sub>S test strip.
16. A system according to any of claims 12 to 15, or a method according to any of claims 1 to 8, wherein said biochemical reagent means include a lysine decarboxylase detection reagent, or arginine decarboxylase and ornithine decarboxylase.
17. A system according to any of claims 12 to 16, or a method according to any of claims 1 to 8, wherein the biochemical reagent means include substrate utilization reactions including as a substrate dulcitol, propylene glycol, and glucuronate.
18. A system according to any of claims 12 to 17, or a method according to any of claims 1 to 8, wherein said biochemical reagent means include methyl umbelliferyl caprylate (MUCAP) and methyl umbelliferyl glucuronide (MUG) (typically in solution or in solid form).

19. A system according to any of claims 12 to 18, or a method according to any of claims 1 to 8, wherein said detection means includes immunological assay means for performing immunological assays on said sample.

20. A system according to any of claims 12 to 19, or a method according to any of claims 1 to 8, wherein said container means for containing said sample and growth medium includes a resilient plastic bag.

21. A methyl umbelliferyl caprylate (MUCAP) reagent for use in diagnostic microbiological analysis, said reagent comprising a solution of methyl umbelliferyl caprylate dissolved in ethyl acetate.

22. A reagent according to claim 21, wherein the concentration of methyl umbelliferyl caprylate in said solution ranges from approximately 0.5% to 10%, preferably approximately 0.1%.

FIGURE 1

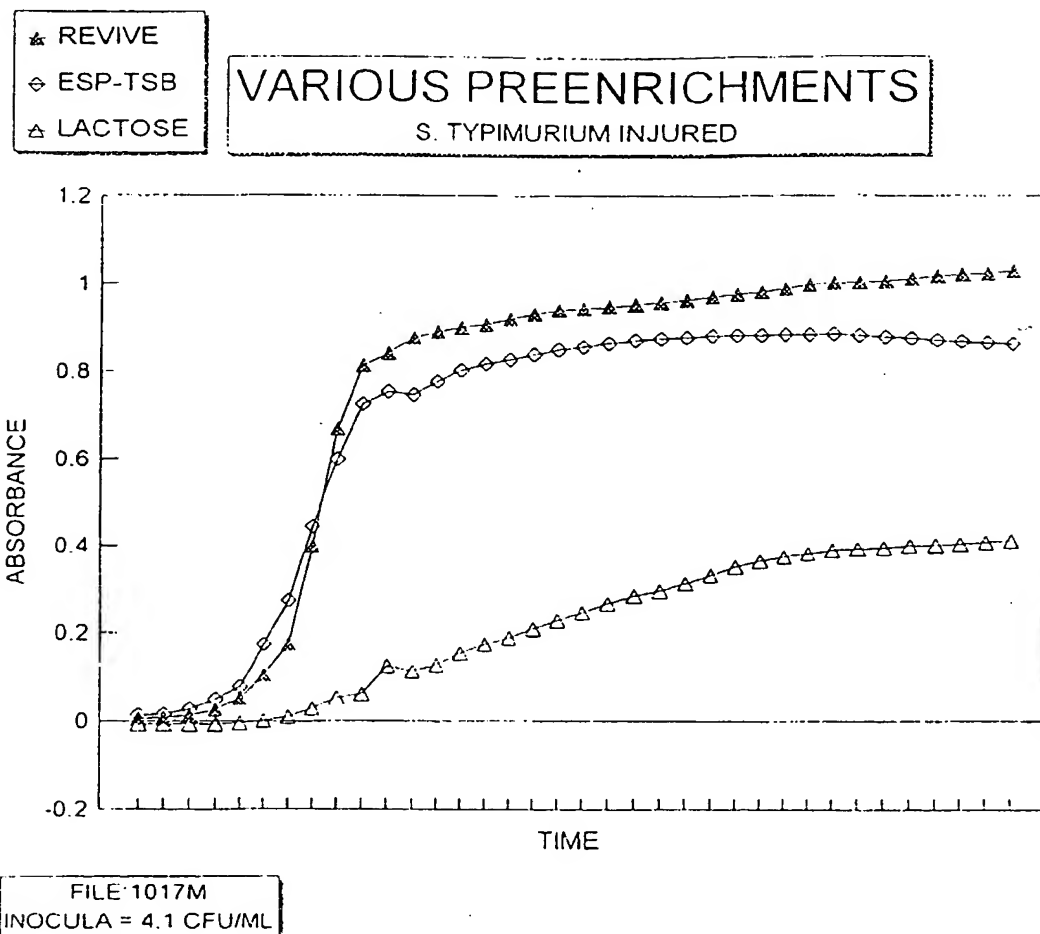


FIGURE 2


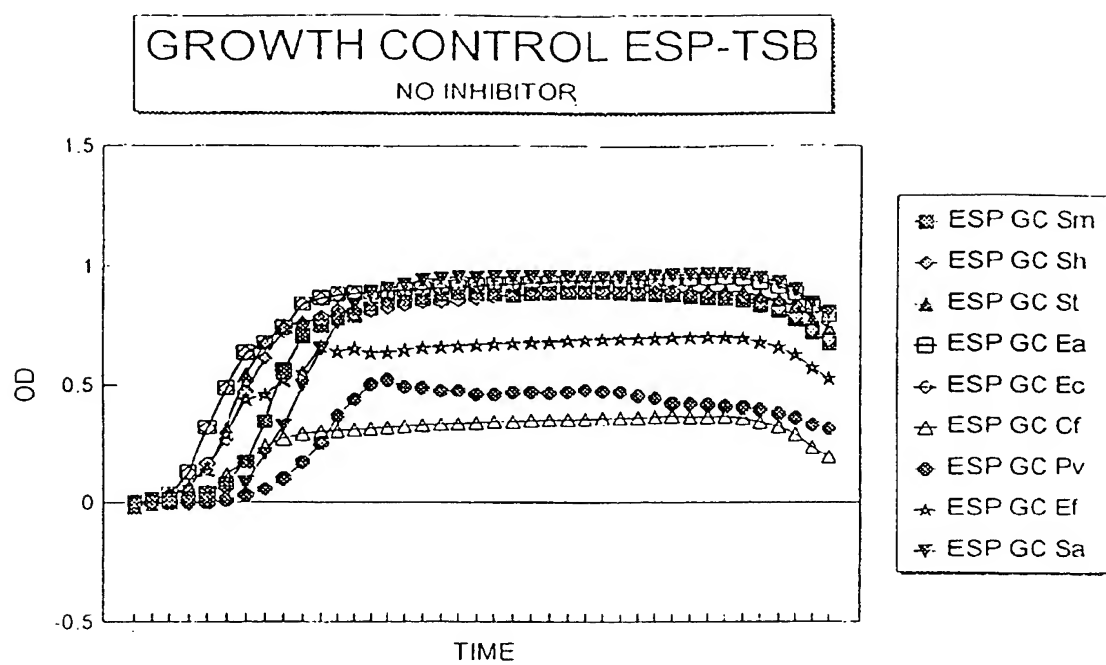


FIGURE 3

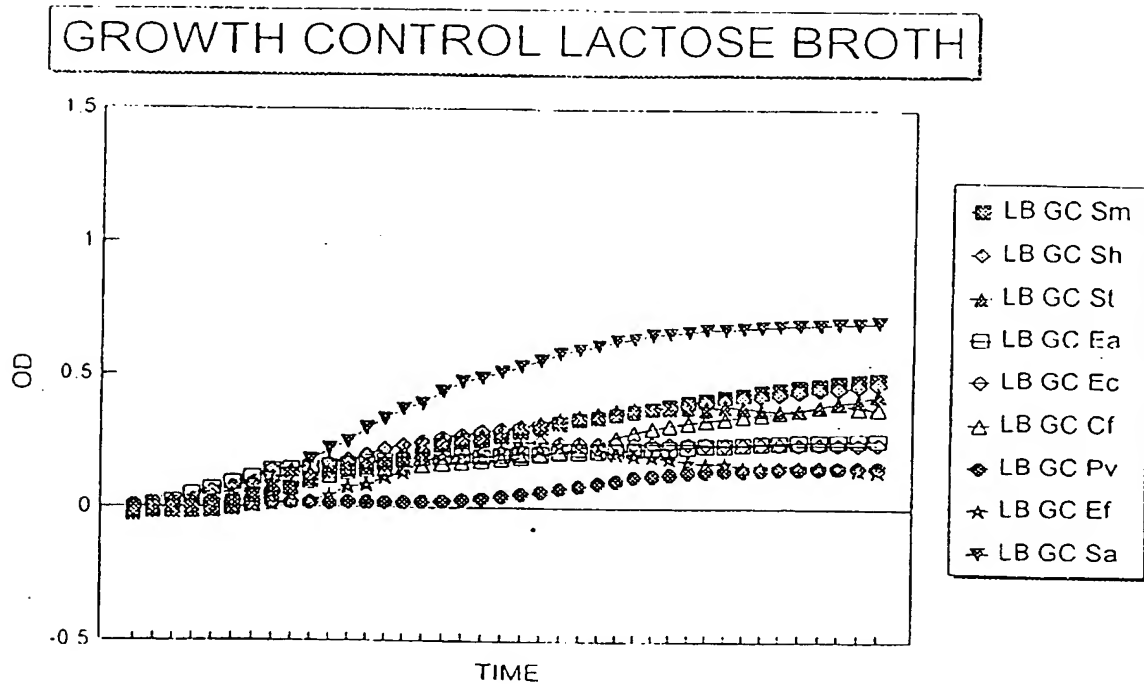

ESP GC= ESP-TSB Aerobic Broth Growth Control

LB GC= Lactose Broth Growth Control

Sm= *S. muenster*, Sh= *S. heidelberg*, St= *S. typhimurium*, Ea= *Enterobacter aerogenes*,

Ec= *E. coli*, Cf= *Citrobacter freundii*, Pv= *Proteus vulgaris*, Ef= *Enterococcus faecalis*

Sa= *Staphylococcus aureus*

FIGURE 4

## ESP WITH BRILLIANT GREEN/DEOXYCHOLATE


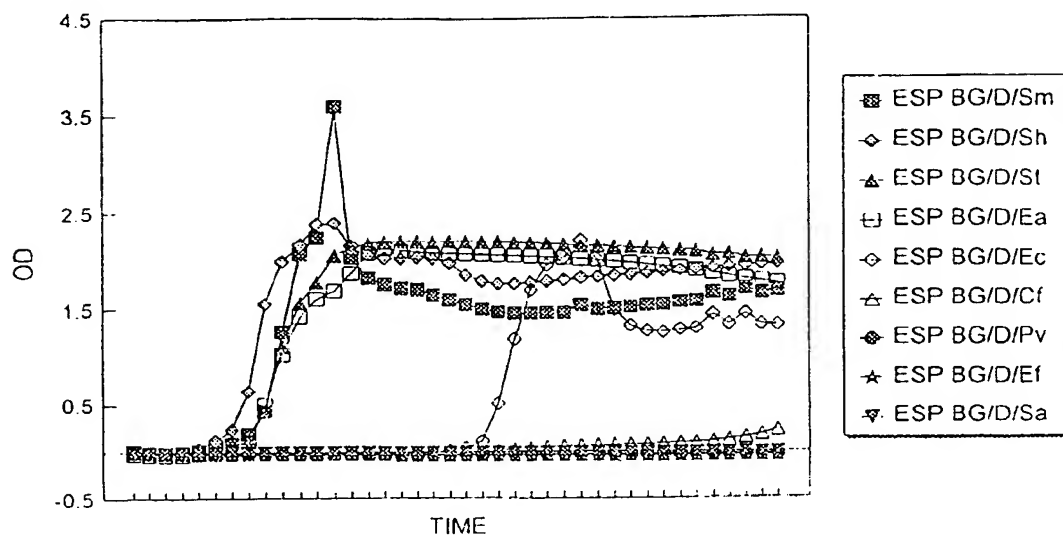


FIGURE 5


## LB WITH BRILLIANT GREEN/DEOXYCHOLATE

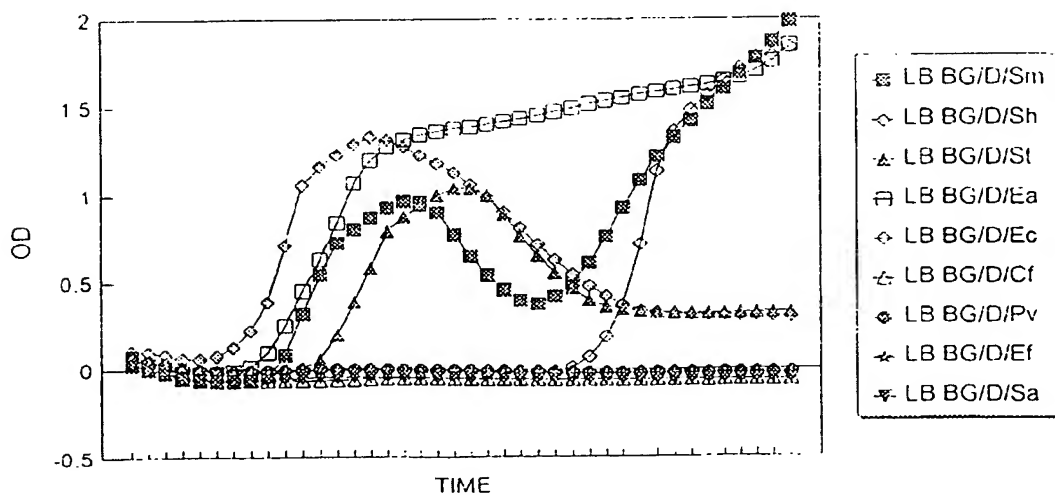

Y-axis has different scale

FIGURE 6

## ESP WITH BRILLIANT GREEN/NOVOBIOCIN

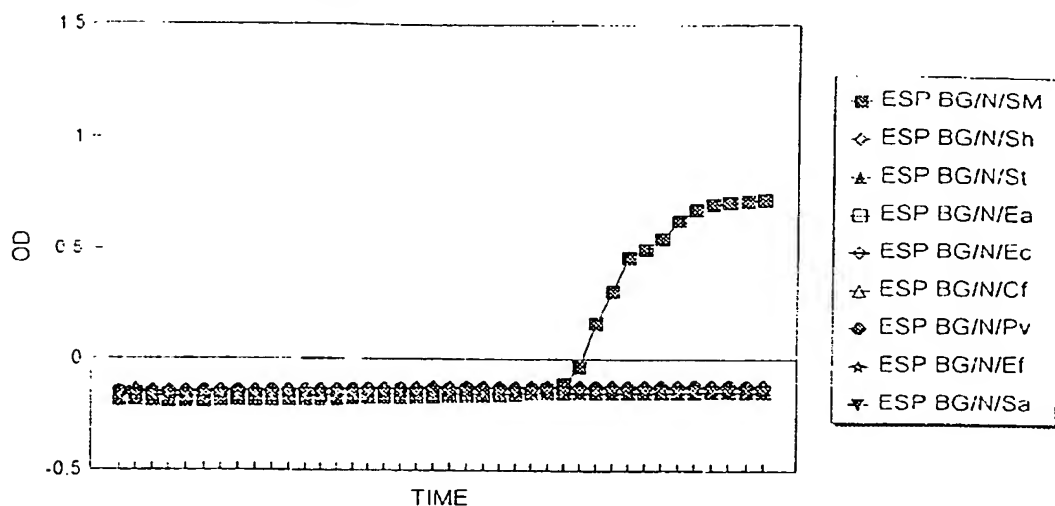

NOTE: LACTOSE BROTH WITH BRILLIANT GREEN/NOVOBIOCIN:  
ALL ORGANISMS WERE INHIBITED

FIGURE 7

## ESP-TSB WITH TERGITOL 4

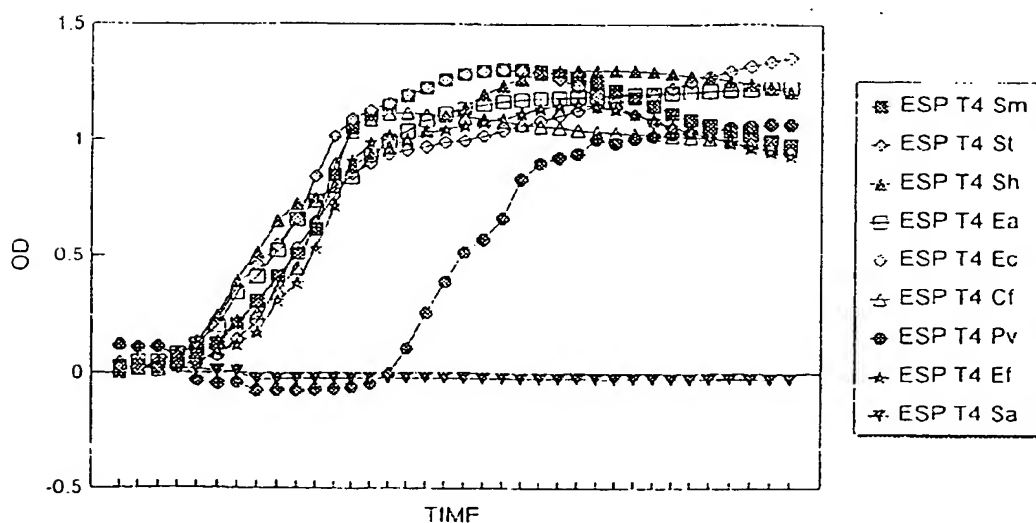


FIGURE 8


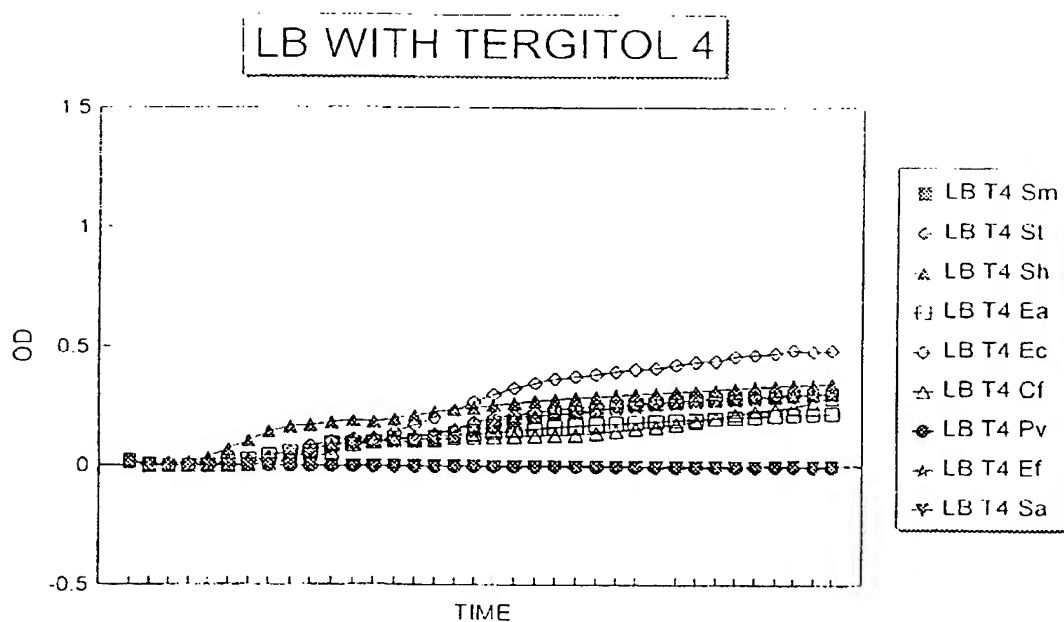




FIGURE 9

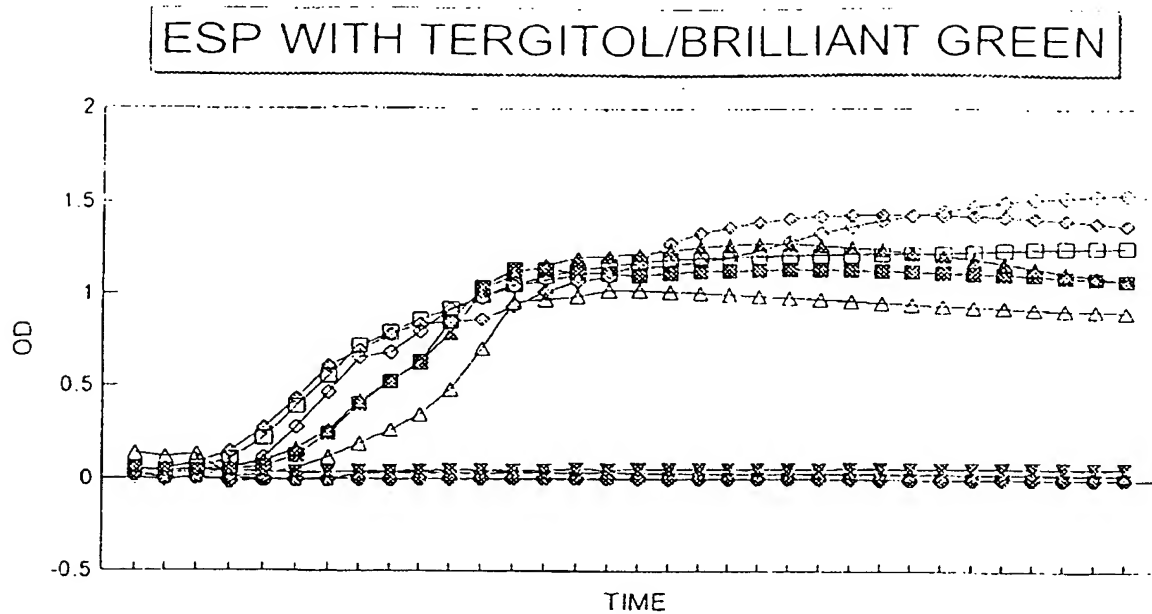

Y-axis has different scale

NOTE: LACTOSE BROTH WITH TERGITOL/BRILLIANT GREEN:  
ALL ORGANISMS WERE INHIBITED

- LB T4 Sm
- ◆ LB T4 Sl
- ★ LB T4 Sh
- ⊖ LB T4 Ea
- ⊕ LB T4 Ec
- △ LB T4 Cf
- ◇ LB T4 Pv
- ☆ LB T4 Ef
- ▽ LB T4 Sa

FIGURE 10


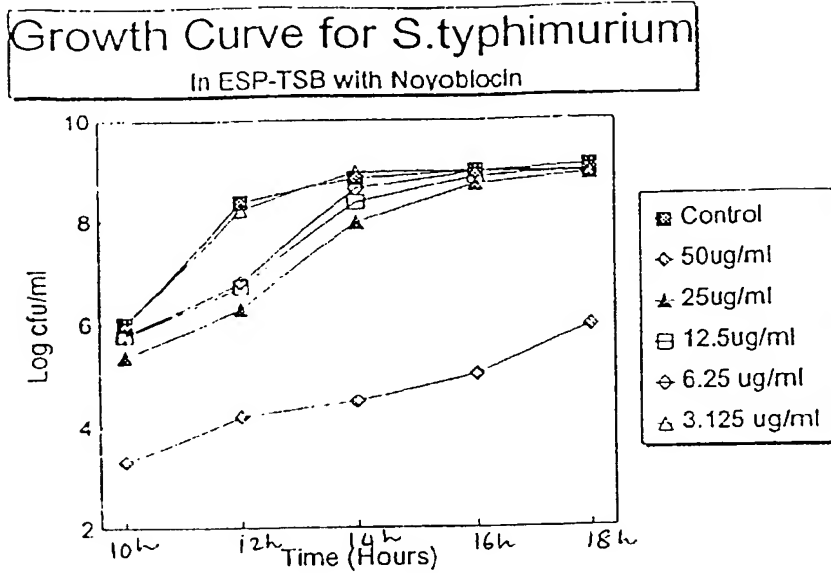


FIGURE 11



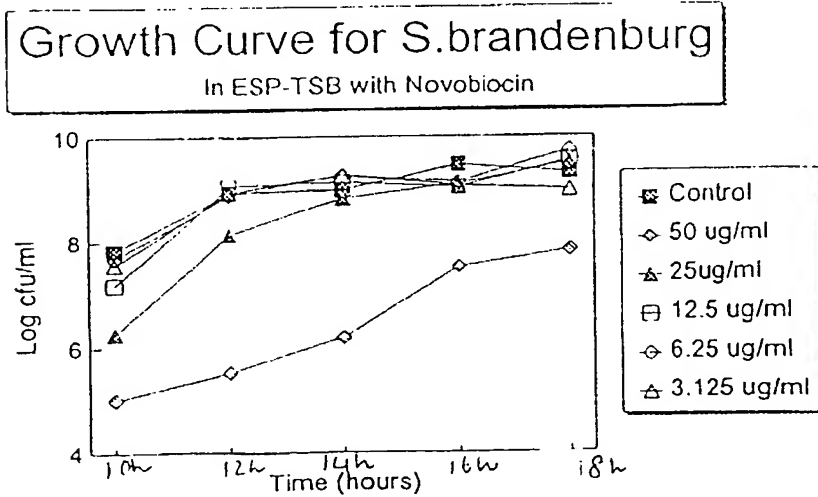

European Patent  
Office

## EUROPEAN SEARCH REPORT

Application Number  
EP 98 30 2542

DOCUMENTS CONSIDERED TO BE RELEVANT			
Category	Citation of document with indication, where appropriate, of relevant passages	Relevant to claim	CLASSIFICATION OF THE APPLICATION (Int.Cl.6)
X	WO 94 28163 A (FOSS ELECTRIC AS :GLENSBJERG MARTIN (DK); HANSEN FLEMMING (DK)) 8 December 1994 • the whole document *	1-13, 15-18	C1201/04
X	WO 92 17607 A (UNIV MARYLAND) 15 October 1992 • abstract *	1,12-14	
X	US 5 447 849 A (TOORA SYED) 5 September 1995 • column 6 - column 8 *	1-4,12	
X	US 4 035 238 A (MEYER MICHAEL C ET AL) 12 July 1977 • abstract *	1,12	
X	US 5 411 867 A (CHANG GEORGE W ET AL) 2 May 1995 • abstract *	1,12	
X	CHEMICAL ABSTRACTS, vol. 117, no. 23, 7 December 1992 Columbus, Ohio, US; abstract no. 229537, MUELLER H. E.: "Afunctional analysis of the components of Preuss tetrathionate broth for the enrichment of Salmonellae" XP002073478 * abstract * & LABORATORIUMSMEDIZIN, vol. 16, no. 4, 1992, pages 149-154, -/--	1-4,12	TECHNICAL FIELDS SEARCHED (Int.Cl.6) C12Q
The present search report has been drawn up for all claims			
Place of search THE HAGUE		Date of completion of the search 3 August 1998	Examiner Moreno, C
<p>CATEGORY OF CITED DOCUMENTS</p> <p>X : particularly relevant if taken alone Y : particularly relevant if combined with another document of the same category A : technological background O : non-written disclosure P : intermediate document</p> <p>T : theory or principle underlying the invention E : earlier patent document, but published on, or after the filing date D : document cited in the application L : document cited for other reasons &amp; : member of the same patent family, corresponding document</p>			

EPO FORM 1503/03-82 (P4C01)



European Patent  
Office

# EUROPEAN SEARCH REPORT

Application Number  
EP 98 30 2542

DOCUMENTS CONSIDERED TO BE RELEVANT			
Category	Citation of document with indication, where appropriate, of relevant passages	Relevant to claim	CLASSIFICATION OF THE APPLICATION (Int.Cl.6)
X	MUÑOZ P. ET AL: "Rapid screening of Salmonella species from stool cultures." AMERICAN JOURNAL OF CLINICAL PATHOLOGY, vol. 100, no. 4, October 1993, pages 404-406, XP002073473 * abstract *	21	
X	RUIZ, J. ET AL: "Modification of the methodology of stool culture for Salmonella detection" JOURNAL OF CLINICAL MICROBIOLOGY, vol. 30, no. 2, February 1992, pages 525-526, XP002073474 * abstract *	21	
X	HUMBERT, F. ET AL: "Rapid identification of Salmonella from poultry meat products by using "Mucap test"." INTERNATIONAL JOURNAL OF FOOD MICROBIOLOGY, vol. 8, no. 1, February 1989, pages 79-84, XP002073475 * abstract *	21	
X	OLSSON, M. ET AL: "Identification of Salmonellae with the 4-methylumbelliferyl caprilate fluorescence test." JOURNAL OF CLINICAL MICROBIOLOGY, vol. 29, no. 11, November 1991, pages 2631-2632, XP002073476 * abstract *	21	
P,X	US 5 650 290 A (GRANT MICHAEL A) 22 July 1997 * the whole document *	1-18	
-/--			
The present search report has been drawn up for all claims			
Place of search <b>THE HAGUE</b>		Date of completion of the search <b>3 August 1998</b>	Examiner <b>Moreno, C</b>
<p><b>CATEGORY OF CITED DOCUMENTS</b></p> <p>X : particularly relevant if taken alone Y : particularly relevant if combined with another document of the same category A : technological background O : non-written disclosure P : intermediate document</p> <p>T : theory or principle underlying the invention E : earlier patent document, but published on, or after the filing date D : document cited in the application L : document cited for other reasons &amp; : member of the same patent family, corresponding document</p>			

EPO FORM 1503 03.02 (P/4C01)



European Patent  
Office

# EUROPEAN SEARCH REPORT

Application Number  
EP 98 30 2542

DOCUMENTS CONSIDERED TO BE RELEVANT			
Category	Citation of document with indication, where appropriate, of relevant passages	Relevant to claim	CLASSIFICATION OF THE APPLICATION (Int.Cl.6)
P, X	WO 97 32034 A (OXOID LTD ; STEPHENS PETER JEREMY (GB)) 4 September 1997 * the whole document *	1, 12	
A	KROMOREDJO, P. & FUJIOKA, R. S.: "Evaluating three simple methods to assess the microbial quality of drinking water in Indonesia." ENVIRONMENTAL TOXICOLOGY AND WATER QUALITY, vol. 6, no. 2, May 1991, pages 259-270, XP002073477 * abstract *	15, 18	
A	DATABASE WPI Section Ch, Week 9248 Derwent Publications Ltd., London, GB; Class B04, AN 92-398887 XP002073479 & WO 92 19764 A (BAXTER DIAGNOSTICS INC) * abstract *	20	
A	US 5 510 243 A (BOYD STEVEN H ET AL) 23 April 1996 * column 3 - column 4 *	1, 12	TECHNICAL FIELDS SEARCHED (Int.Cl.6)
A	US 4 279 995 A (WOODS LEODIS V ET AL) 21 July 1981 * abstract *	1, 12	
A	US 5 449 612 A (LEPARGNEUR JEAN-PIERRE ET AL) 12 September 1995 * abstract *	1, 12	
A	WO 96 06182 A (MINNESOTA MINING & MFG) 29 February 1996 * abstract *	1, 12	
The present search report has been drawn up for all claims			
Place of search THE HAGUE		Date of completion of the search 3 August 1998	Examiner Moreno, C
<p>CATEGORY OF CITED DOCUMENTS</p> <p>X : particularly relevant if taken alone Y : particularly relevant if combined with another document of the same category A : technological background O : non-written disclosure P : intermediate document</p> <p>T : theory or principle underlying the invention E : earlier patent document, but published on, or after the filing date D : documents cited in the application L : document cited for other reasons &amp; : member of the same patent family, corresponding document</p>			

EPO FORM 1503 03 82 (P4/C01)

**THIS PAGE BLANK (USPTO)**

**This Page is Inserted by IFW Indexing and Scanning  
Operations and is not part of the Official Record**

**BEST AVAILABLE IMAGES**

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images include but are not limited to the items checked:

☐ **BLACK BORDERS**

☐ **IMAGE CUT OFF AT TOP, BOTTOM OR SIDES**

☐ **FADED TEXT OR DRAWING**

☒ **BLURRED OR ILLEGIBLE TEXT OR DRAWING**

☐ **SKEWED/SLANTED IMAGES**

☐ **COLOR OR BLACK AND WHITE PHOTOGRAPHS**

☐ **GRAY SCALE DOCUMENTS**

☒ **LINES OR MARKS ON ORIGINAL DOCUMENT**

☐ **REFERENCE(S) OR EXHIBIT(S) SUBMITTED ARE POOR QUALITY**

☐ **OTHER:** \_\_\_\_\_

**IMAGES ARE BEST AVAILABLE COPY.**

**As rescanning these documents will not correct the image problems checked, please do not report these problems to the IFW Image Problem Mailbox.**

**THIS PAGE BLANK (USPTO)**